




*Making Life easier*



*Solutions to Daily Living*

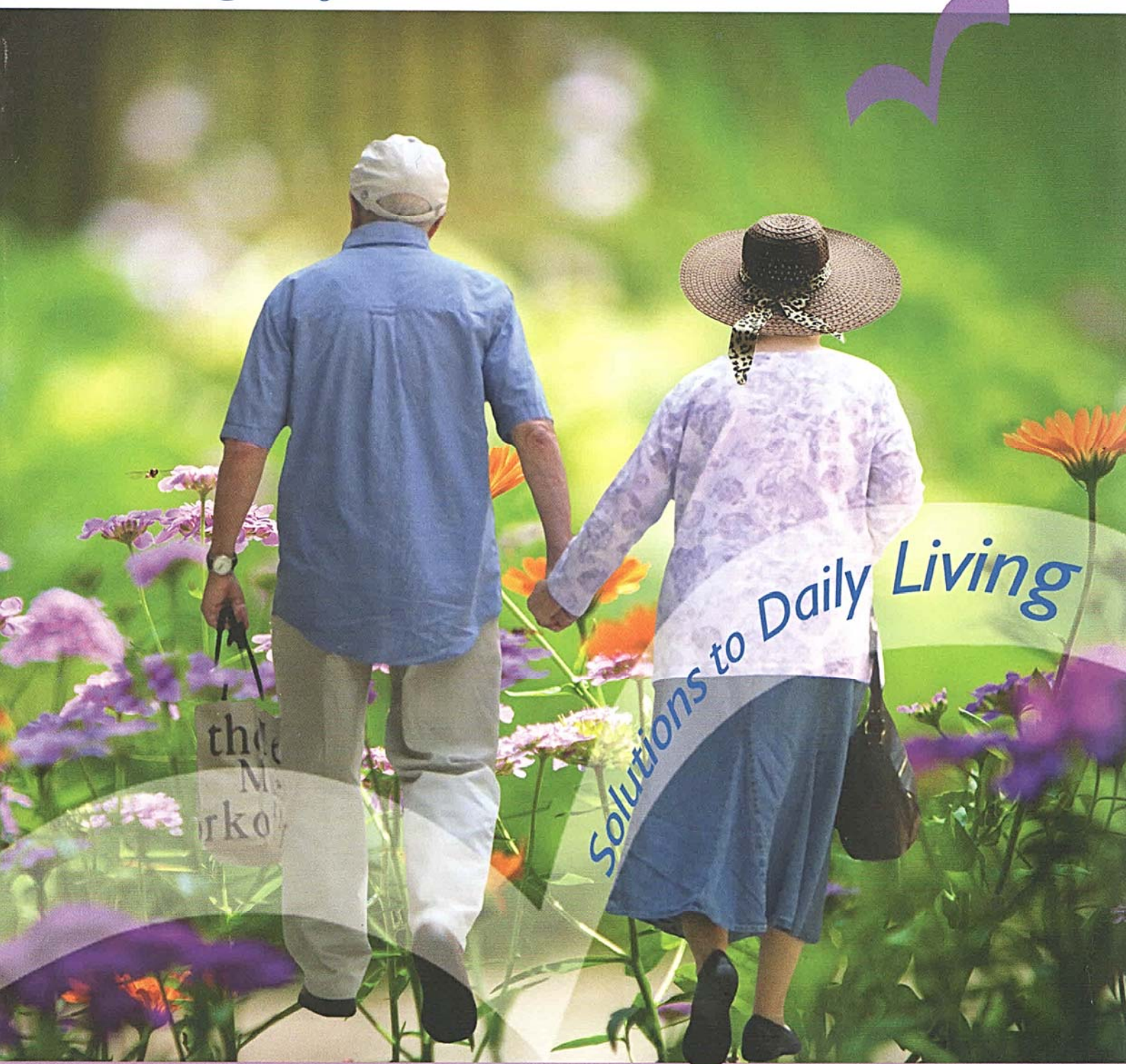

+961 3 711112 / +961 9 224004 [www.homemobility.info](http://www.homemobility.info)  
[info@homemobility.info](mailto:info@homemobility.info)

## Introduction

### INDEX

Seating	Page 1 - 6
Riser & Recliner	Page 2 - 5
Recliner & Cushion	Page 6
Bedroom	Page 7 - 13
Bed	Page 8 - 9
Mattress	Page 10 - 11
Accessories	Page 12
Bathroom	Page 13 - 18
Bathing	Page 14 - 15
Toileting	Page 16
Showering	Page 17
Special Design	Page 18
Grab Bar	Page 18
Access	Page 19 - 28
Hoist	Page 20 - 21
Ramp	Page 22
Stairlift	Page 23 - 25
Platform Lifts	Page 26
Stair Climber	Page 26
Elevator	Page 27
Regain your Freedom	Page 28
Mobility	Page 29 - 36
Walking Aid	Page 30
Wheelchair	Page 31 - 35
Cushions and Backs	Page 36
Easy Stand	Page 36

At **Home Mobility**, we provide Solutions to everyday mobility problems in and outside home. Our goal is to improve your Quality of Life.

We start with a Free Assessment & we make sure to give you a solution tailored to your specific needs. Our Qualified Team is professionally trained to give you the right information.

Our products are designed to promote an independent and involved life style: they are carefully selected from reputed brands & meet the Highest Safety Standards.

We do lot of research to provide you with innovative programs:



We developed a flexible Rental Program for your short term need. All rented items are cleaned, disinfected and checked after each usage.



We can tailor some of our products to your specific needs.

Purchased items are delivered free of charge, anywhere in Lebanon; in addition, we provide Warranty, After-Sale Services & Consultancy.

Contact us at: 961-(0)3-711 112 / 961-(0)9-224 004  
[info@homemobility.info](mailto:info@homemobility.info) / [www.homemobility.info](http://www.homemobility.info)  
 P.O.Box: 428 Zouk Mikael - Lebanon

or visit us in our showroom:  
 Villa Dahan, T2 Street, Old Souk East Entrance,  
 near Municipality, Zouk Mikael, Lebanon.

## Making Life Easier, that's what we do!

Home Mobility Team

We are proud of OUR BRANDS



مرحباً!  
 اسمي أيمن وهدفي سهلكن  
 عيشتمكم.



*Seating*



Home Mobility Lebanon



Lift Chair Lebanon

[www.homemobility.info](http://www.homemobility.info)

# Seating

## Riser & Recliner



### Riser & Recliner®



Sleep



Read



Sit



Stand

**Riser & Recliners** help you stand up and recline to sleep-ing position effortlessly, at the push of a button. In addition, they are luxurious & comfortable chairs, Hand Crafted using the highest quality components to give you maxi-mum service with minimum maintenance.

Stand up safely  
at the push of a button!

Only genuine Golden Riser & Recliner® offers the Golden Warranty



Hand-blown, Overstuffed,  
Conjugated Microfiber Dacron Fill  
Super Supportive

Quick Release Hand Control  
Easy Disconnect


Widest Fabric Selection  
Fits in any Home Décor

5 Heavy Gauge Springs  
Long Lasting for Incredible Comfort

10 Gauge Steel Scissor Mechanism  
Strong & Long Lasting



Battery backup  
ETL approved

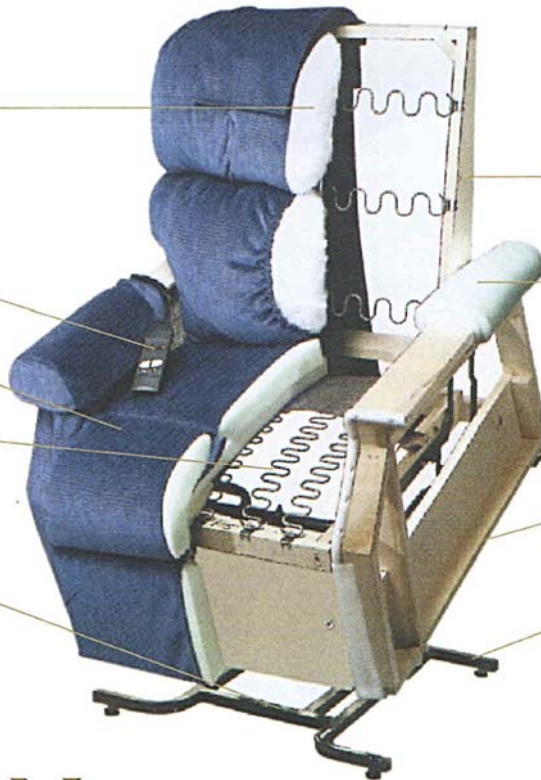

Northern Maple Wood  
-Double-doweled at Stress Points  
- Reinforced Miter Blocks  
-Stapled & Glued  
Strongest in the industry

2" of Block Foam  
Extra Comfort

Rear Landing Legs  
Extra Stability

Okin Motor  
Most Reliable and Quiet

Adjustable Anti-Skid/Scrape Feet  
Stability & Support

# GOLDEN

technologies

Life time warranty on Lift & Chair Frames and Recliner Mechanism.  
One year warranty on defective parts, excluding misuse, wear & tear  
or improper use.

2 | ☎ 03 - 711112 | 09 - 224004



عندك صعوبة بالنهوض من الكرسي؟  
الحل هو هلكرسي:  
بتساعد على النهوض وبتتحول إلى سرير.  
كبسة زر وإرتاح مثل ما بذك!



## Golden Comforter Series



### Made to Measure

Seat Dim. wxdxh:  
O/A Dim. wxdxh:  
Distance from wall:  
Max. user weight:  
User height:



### Ref. PR501/L

58 x 56 x 52 cm  
88 x 83 x 112 cm  
46 cm  
170 Kg  
178 - 187 cm



### Ref. PR501/M

53 x 53 x 50 cm  
83 x 83 x 107 cm  
43 cm  
170 Kg  
162 - 177 cm



### Ref. PR501/P

48 x 47 x 46 cm  
77 x 83 x 99 cm  
43 cm  
136 Kg  
< 152 cm

**GOLDEN**  
technologies



## Comforter

Tired of trying to fit in a Riser & Recliner that's sized for someone else? Look to the Comforter Series to provide your optimum fit. The Comforter's plush design provides exceptional lumbar support and a full chaise pad for total body comfort.



Admiral



Palamino



Evergreen



Cabernet

## Golden Value Series



Rental Available



Vinyl



Autumn



Cornflour



Chestnut



Spring



Rosewood



**New**

## Capri 3

This Riser & Recliner chair has a classic colonial button back design enhanced by a slightly firmer lumbar support area. The arms offer a wider, more rounded design for extra comfort and stability.

### Ref. PR350-3

Seat Dim.: 51 x 48 x 46 cm - O/A Dim.: 79 x 83 x 104 cm  
Distance from wall: 43 cm - Max. user weight: 147 Kg

## Monarch Plus

The Monarch is a distinctive Riser & Recliner that combines a generously stuffed seam back with a soft supportive full chaise seat. Combines ample room to relax and a tasteful design.

### Ref. PR359-M

Seat Dim.: 53 x 55 x 50 cm - O/A Dim.: 81 x 83 x 105 cm  
Distance from wall: 46 cm. Max. user weight: 170 Kg

## Golden Traditional Series



Ref. PR747

Seat Dim.: 53 x 55 x 50 cm  
O/A Dim.: 83 x 83 x 112 cm  
Distance from wall: 46 cm  
Max. user weight: 170 Kg



Ref. PR643

Seat Dim.: 55 x 53 x 48 cm  
O/A Dim.: 75 x 83 x 104 cm  
Distance from wall: 46 cm  
Max. user weight: 170 Kg



Galaxy



Burgundy



Pine



Black



Cream



Island



Sterling



Caramel



Indigo

### Williamsburg

Wrap yourself in comfort! The Williamsburg is a Riser & Recliner that features a soft colonial design with traditional wings. This fashionable button-back design provides comfort and proper support for much needed relaxation.

### Pioneer

Easy to grasp wooden arms will provide you with more stability while getting in or out of this chair. The seat and back cushions will provide sink in softness for levels of comfort all its own. With its youthful appeal, the Pioneer is the ideal complement for any home or office.

## Golden Signature Series



### Space Saver

The Space Saver is a unique Riser & Recliner in that it only requires 15 cm of space from the wall to recline. The overstuffed back and arms will invite you to relax for hours of satisfying comfort. Include heat and vibration massage.

Ref. PR906

Seat Dim.: 51 x 56 x 56 cm  
O/A Dim.: 90 x 83 x 109 cm  
Distance from wall: 15 cm  
Max. user weight: 113 Kg  
User height: > 180 cm



### Regal

The Regal is a feature rich Riser & Recliner designed for people that like to keep things within hands reach. The dual flip open arms conceal extra storage space, a cup holder and hideaway tray. The luxurious three pillow waterfall back can be adjusted, providing more comfort options.

Ref. PR751

Seat Dim.: 56 x 55 x 51 cm  
O/A Dim.: 91 x 83 x 111 cm  
Distance from wall: 46 cm  
Max. user weight: 170 Kg

### FABRIC



Saddle



Havana



Vino



Orion

## Golden Maxi Comfort Series



Golden Technologies has perfected comfort. With the development of our exclusive, patented MaxiComfort® positioning technology, Golden delivers the most advanced dual motor Riser and Recliner chair available anywhere in the world!

experience the  
*Health Benefits*

*Rejuvenate*  
with a great night's sleep

*Relieve Stress*  
with the pressure-free luxury of zero-gravity!

*Relax & Watch TV*  
in comfort




### Relaxer

Four-pillow waterfall back, adjustable lumbar support & plush rounded armrests makes you feel like you're resting on a cloud.

Ref. PR756

Seat Dim.: 51 x 53 x 48 cm  
O/A Dim.: 93 x 91 x 111 cm  
Distance from wall: 64 cm  
Max. user weight: 170 Kg

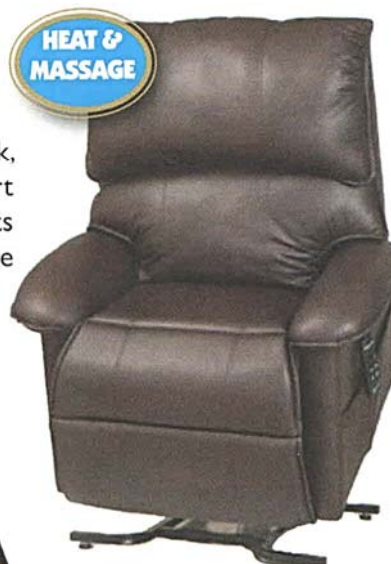


### Windsor

Two-pillow backrest design with contemporary, piecrust seams enhance the armrests and backrest design. Its casual style blends with any décor.

Ref. PR516

Seat Dim.: 53 x 51 x 53 cm  
O/A Dim.: 99 x 86 x 114 cm  
Distance from wall: 71 cm.  
Max. user weight: 170 Kg

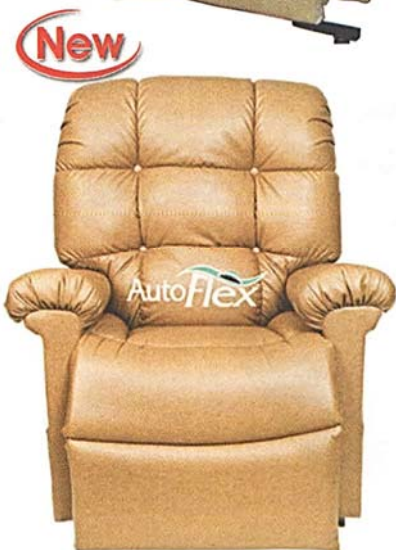

### Cloud

The most comfortable Riser & Recliner Golden have ever built! An overstuffed biscuit style back and a unique seat that cradles you in comfort. AutoFlex power lumbar system gently stretches the spine & muscles in the lower back.

Ref. PR510

Seat Dim.: 51 x 55 x 51 cm  
O/A Dim.: 97 x 90 x 114 cm  
Distance from wall: 66 cm.  
Max. user weight: 170 Kg



**AutoDrive**  
CUSTOMIZED COMFORT



### Cirrus Large

Features Golden's pocketed coil spring seat, surrounded by medium density foam, molds to your body for super soft, yet supportive comfort ... the seat stays with you as you shift your weight from side to side.

Ref. PR508L

Seat Dim.: 56 x 53 x 55 cm  
O/A Dim.: 90 x 88 x 117 cm  
Distance from wall: 71 cm.  
Max. user weight: 170 Kg

**BRISA®**




Buckskin Coffee Bean Sepia Black Onyx

**New**

**GOLDEN**  
technologies

BRISA Ultra-leather is an incredible soft & luxurious faux leather polyurethane fabric. Its unique ventilation system enables it to breathe, keeping you cool, no matter how long you use the chair. BRISA FRESCO is enhanced with silver ion antimicrobial agents.

## Recliner

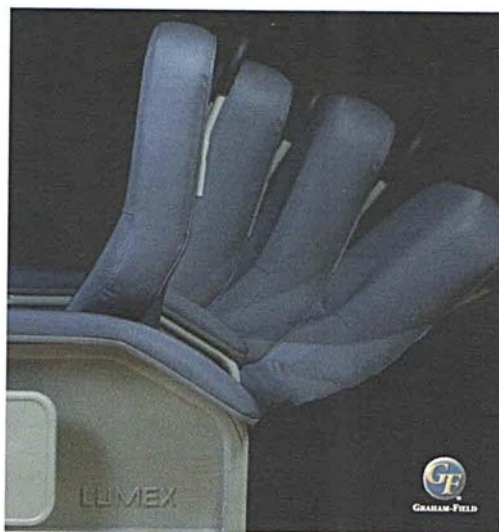

### Preferred Care Recliner with drop arms


Infinite recliner back from 95° to 145° and a unique suspension system cradle the user from head to heel for increased comfort and pressure relief. Self-storing, height adjustable drop arms enables a safer & easier lateral transfer. Directional wheel lock on right rear caster reduces "crabbing" of recliner. Quick-release upholstery allows for easy cleaning, removal and replacements. Vinyl upholstered back is sulfide stain, mildew and oil resistant.

Ref. 565DG

O/A Dim: 71 x 91 x 122 cm

Max. user weight: 113 Kg

## Cushion

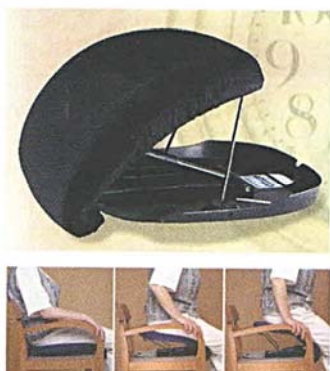

### Uplift Seat Assist

Stand up with this self-powered cushion. Take it anywhere you go.

Ref. UL 100 / UL 300

Dim.: 43 x 46 x 7 - 24 cm (WxDxH)

Max. user weight: 35 - 105 / 90 - 160 Kg

## Raising Blocs



COOPERS

### Raising Blocs

Raise your bed or arm chair; it will help you stand up.

Ref. 10720 / 10711

Height: 5 - 10 cm

Max. user weight: 190 Kg



### Permalux Cushion

Economical cushion, fiber-filled with waterproof, fire retardant cover.

Ref. PL024

O/A Dim: 43 x 43 x 8 cm



### Permapad Cushion

CMHR modular cut foam seat with bi-stretch, zipped, waterproof, fire retardant cover.

Ref. PL077 / PL078

O/A Dim: 43 x 43 x 10 cm

51 x 51 x 10 cm

Max. user weight: 108 Kg



### Memaflex Cushion

Visco-elastic foam seat conforms to the body with bi-stretch, zipped, waterproof fire retardant cover

Ref. PL106

O/A Dim: 43 x 43 x 7.5 cm

Max. user weight: 127 Kg



# Bedroom



Home Mobility Lebanon

[www.homemobility.info](http://www.homemobility.info)

## Beds



Rental  
Available

### Patriot Semi-Electric Bed

The robust grid sleep surface, made in USA, provides better mattress support than traditional spring decks. Head & foot/knee-break sections can be raised or lowered simultaneously or independently with a hand control. The bed height can be adjusted manually with a crank. Okin motor provides a quiet, smooth operation, with a 9v battery back-up. Bed includes 4 castors (2 lockable) & fiberboard decorative walnut bed ends.

Ref. US0218

Elevation: 60° head, 30° knee

O/A Dim: 91 x 220 x 33 – 58.4 cm

Max. user weight: 158 Kg (204 Kg w/accessories)



Rail Slot



Locking  
Pins



Rental  
Available

### Patriot Full-Electric Bed

This bed offers the same features as the Patriot Semi-Electric Bed with the additional feature of electric height adjustment. It makes it easier for caregivers to raise and lower the bed for patient care and linen changes. The electric height adjustment also enables the patient to easily get in and out of the bed with a maximum comfort.




Ref. US0458

Elevation: 60° head, 30° knee

O/A Dim: 91 x 220 x 33 – 58.4 cm

Max. user weight: 158 Kg (204 Kg w/accessories)

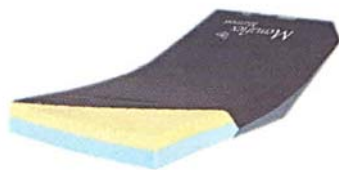

### Mattress

Select Your Mattress on Page 11.



### Silver Solutions

Powerful antimicrobial technology that suppresses the growth of bacteria.



### Liberty Bed Rails

Designed to offer a safe and secure bed rail system.



Special  
Order

### Forza Bariatric Bed

Heavy duty full electric bed with 4-section steel mesh bed surface, available in various width and weight configuration. Easy transport despite size or weight. Optional: bed to chair position

Ref. Forza

O/A Dim: 111 / 151 / 35w x 200L x 40 – 80h cm

Max. user weight: 250 / 350 / 500 Kg



سرير المستشفى بالبيت:  
يساعد على النهوض وبسهولة الاهتمام  
بالمريض.  
تمتع بنوم هنيء!



## beds



Long term  
care

Special  
Order



Central lock



Electrical Safety Plug



High Quality Ratchet



Pendant with Lock

### Allegra Full Electric Bed

Full electric bed with homely & elegant design includes wooden headboard & footboard with integrated security side rails. Very sturdy four section steel mesh on metallic feet with integrated double castors and central lock. Lockable pendant operates the back rest, leg rest with knee-break, bed height, Trendelenburg & Reverse-Trendelenburg functions. Double-drive Okin motors ensures ultra-silent bed operation, with 9V battery back-up.

Ref. 2EN4083

Elevation: 70° head, 45° knee

Trandelenburg: + / - 11°

O/A Dim: 105 x 208 x 40 – 80 cm

Max. user weight: 155 Kg (170 Kg w/accessories)

Made in  
Germany



Fowler



Vascular



Lo



Hi



Trendelenburg



Cardiac Chair



Long term  
care

Special  
Order

Ref. 2M4081

Elevation: 70° head, 20° knee

Trandelenburg: + / - 11°


O/A Dim: 103 x 211 x 40 – 80 cm

Max. user weight: 135 Kg (170 Kg w/accessories)

mühle  
müller pflegebetten

### Medial

Full electric bed with wooden headboard & footboard and integrated security side rails. Two Dewert motors electrically adjust the back rest, legrest with knee brake, bed height, Trendelenburg & reverse-Trendelenburg functions. Handset has a magnetic swipe tab to lock functions. 4 lockable castors.

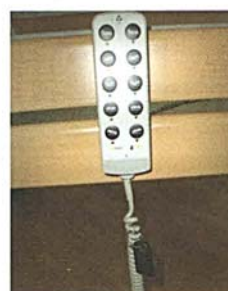

Pendant with  
Magnetic Lock



Transport

## Mattress

### Skin / Wound Assessment

Patient Name: \_\_\_\_\_  
Diagnosis: \_\_\_\_\_ Predisposing Diseases: \_\_\_\_\_ Weight: \_\_\_\_\_  
Date \_\_\_\_\_ Date \_\_\_\_\_

Doctor Name: \_\_\_\_\_  
Instructions: Assess upon admission, every Wednesday thereafter & when condition changes. Inspect skin at least once per day

#### General Physical Condition (Health Problem)

Good	- Usually free of major health problems	0	<input type="text"/>	<input type="text"/>
Fair	- Major health problems but stable	1	<input type="text"/>	<input type="text"/>
Poor	- Chronic / Serious health problems, not stable	2	<input type="text"/>	<input type="text"/>

#### Level of Consciousness to commands

Alert	- Respond readily	0	<input type="text"/>	<input type="text"/>
Lethargic	- Slow to respond	1	<input type="text"/>	<input type="text"/>
Semi-comatose	- Responds only to painful stimuli	2	<input type="text"/>	<input type="text"/>
Comatose	- No response to stimuli	3	<input type="text"/>	<input type="text"/>

#### Activity

Ambulant without assistance		0	<input type="text"/>	<input type="text"/>
Ambulant with Assist of Person/Walker/Prosthesis		2	<input type="text"/>	<input type="text"/>
Chairfast	- Cannot ambulate/Restricted to wheelchair	4	<input type="text"/>	<input type="text"/>
Bedfast	- Unable to sit/Bed-Bound	6	<input type="text"/>	<input type="text"/>

#### Mobility (Extremities)

Full active range	- Able to move extremities at will	0	<input type="text"/>	<input type="text"/>
Restricted Movement	- Moves with assistance	2	<input type="text"/>	<input type="text"/>
Moves only with assistance		4	<input type="text"/>	<input type="text"/>
Immobile		6	<input type="text"/>	<input type="text"/>

#### Incontinence (Bowel and/or Bladder)

None		0	<input type="text"/>	<input type="text"/>
Occasional	- Less than 2 per 24 hrs.	2	<input type="text"/>	<input type="text"/>
Usually	- More than 2 per 24 hrs.	4	<input type="text"/>	<input type="text"/>
Total	- No control	6	<input type="text"/>	<input type="text"/>

#### Nutrition (for age and size)

Good	- Eats/Drinks 3/4 - all of meal or Albumin $\geq 3.5\text{gm/dl}$	0	<input type="text"/>	<input type="text"/>
Fair	- Eats/Drinks 1/2 - 3/4 - all of meal or on TPN or Albumin 3.0-3.4gm/dl	1	<input type="text"/>	<input type="text"/>
Poor	- Unable/Refuses < 1/2 of meal or Albumin 3.0gm/dl	2	<input type="text"/>	<input type="text"/>

#### Skin/Tissue Status

Good	Well nourished/Skin intact	0	<input type="text"/>	<input type="text"/>
Fair	Poorly nourished/Skin intact	4	<input type="text"/>	<input type="text"/>
Poor	Skin not intact	6	<input type="text"/>	<input type="text"/>

IF THE TOTAL SCORE IS 8+ INITIATE SKIN/WOUND PROTOCOL

Low Risk: 1 - 7 Medium Risk: 8 - 22 High Risk: 23 - 31

Score

Carer Name: \_\_\_\_\_


Signature \_\_\_\_\_

#### Describe Skin/Wound Condition & Stage: \_\_\_\_\_

Grade: \_\_\_\_\_

Use Diagrams to note Locations:

- 0 - Skin intact, free of redness
- I - Reddened area
- II - Blister, skin break
- III - Skin break exposing subcutaneous skin
- IV - Skin break exposing muscle, bone or tendons

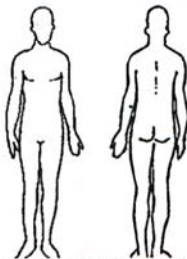

#### Note:

This Information is the property of Home Mobility s.a.r.l. & may not be reprinted without prior permission.  
This form is used at St. Mary Hospital - New York

#### Skin/Wound

#### Protocol

#### Intervention

- Pressure Reduction
- Pressure Relief
- Incontinence Care
- Nutrition Consult
- Moisturize Skin
- Positioning
- Wound Care




### How to select the right mattress

1- Complete assessment form




2- Use the score to select the mattress

3- Follow skin-wound protocol

### Causes of pressure sores

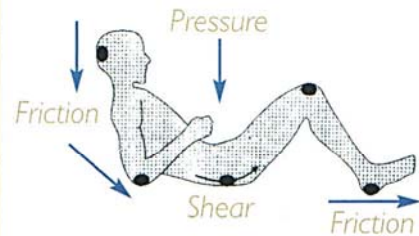

### Positioning

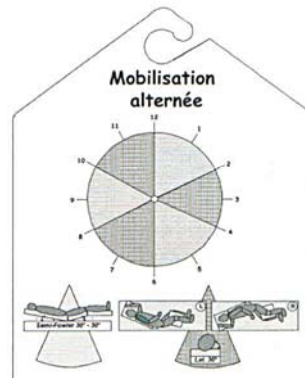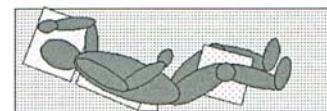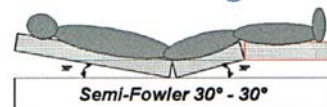

35

30

High risk

25

20

Medium risk

15

10

5

Low risk

#### Pressure Reduction

Pressure Reduction mattress / overlay are made of special foam type that takes the shape of body, thus reducing the pressure. Designed for people having some mobility.

#### Pressure Relief

Pressure Relief mattress/ overlay are made of air cells that gently inflates and deflates in alternance, thus reducing the duration at the points of pressure. Designed for people with reduced mobility.

العقر بصر مشكلة كبيرة اذا ما  
تتصرف بسرعة. كل شخص يمضي معظم  
وقته بالسرير لازمه حماية. اتصل فينا  
لنعطيك المعلومات اللازمة!



## Mattress

### Pressure Reduction

#### Air Layer Profiling Mattress

Same as Memaflex with a top air layer reduces heat & vapour transfer improving the micro-climate between the user and mattress.


Ref. PL115  
O/A Dim: 88 x 198 x 16 cm  
Max. user weight: 248 Kg



#### Memaflex Profiling Mattress / Overlay

4 cm of contoured visco-elastic foam surface allows greater body immersion thus providing greater pressure re-distribution. 12 cm of CMHR base provides further support to improve patient comfort.

Ref. PL111 Mattress / PL110 Overlay  
O/A Dim: 88 x 198 x 16 / 9 cm  
Max. user weight: 248 Kg / 114 Kg

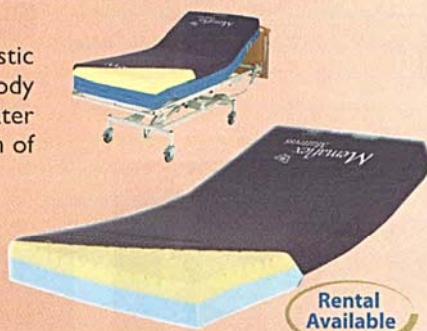

#### Permaflex Profiling Mattress

High degree of comfort & performance. CMHR foam conforms to all body shapes, maximizes weight distribution, recovers excellently & increase air circulation via special cut of mattress.

Ref. PL100  
O/A Dim: 90 x 198 x 16 cm  
Max. user weight: 165 Kg



#### Permapad Profiling Overlay

Cost effective one-piece CMHR foam.

Ref. PL082  
O/A Dim: 88 x 198 x 9 cm  
Max. user weight: 108 Kg



#### Community Profiling Mattress

CMHR foam conforms to all body shapes, maximizes weight distribution, recovers excellently & increases air circulation via special cut of mattress.

Ref. PL060  
O/A Dim: 88 x 198 x 12.5 cm  
Max. user weight: 140 Kg



#### Permalux Profiling Overlay

Fiber filled overlay with stitched cover. Recommended for frail and skinny people

Ref. PL020  
O/A Dim: 90 x 198 x 13 cm  
Max. user weight: 100 Kg




### Pressure Relief

#### Phase III Profiling Air Mattress

Auto Pressure Adjustment to patient movement, rapid inflation, Auto Shut-off valves to prevent mattress deflation during power failure, CPR; quiet pump.

Ref. PA7900  
O/A Dim: 88 x 190 x 20 cm /  
Pump: 12.5 x 35 x 26 cm  
Max. user weight: 29 - 248 Kg

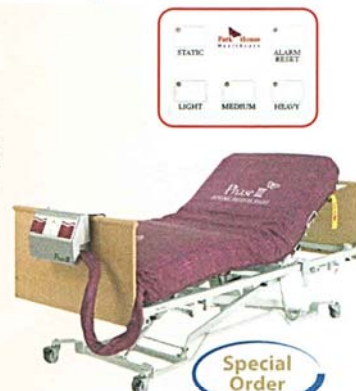

#### Elite / Elite Comfort Profiling Air Mattress

Alternating Air Therapy, 17 individual cells with CMHR foam base prevents bottoming out. CPR, Quiet pump. Elite Comfort has a top Visco layer to enhance patient comfort.

Ref. PA5550/PA5550-79  
O/A Dim: 88 x 200 x 17.5 cm  
Pump: 29 x 20 x 11.7 cm  
Max. user weight: 190 Kg



#### Eclipse II Profiling Overlay

Alternating Air Therapy; 16 individual cells; CPR, quiet pump with audio visual alarm. Include alternating seat.

Ref. PA5500  
O/A Dim: 88 x 185 x 12 cm  
Pump: 20 x 11.7 x 29 cm  
Max. user weight: 28.6 - 139 Kg



#### Prelude Profiling Overlay

Alternating Air Therapy; 18 individual cells; CPR. Quiet pump with visual alarm. Optional alternating seat.

Ref. PA3500  
O/A Dim: 88 x 194 x 10 cm  
Pump: 12.5 x 8.5 x 20 cm  
Max. user weight: 28.6 - 100 Kg



The polyurethane fire-retardant cover, designed to prevent friction & shear, is waterproof, multi stretch, vapour-permeable, impregnated with Acti-fresh biocide and complies with infection control protocol. It can be removed and laundered up to 80° C.



**Park House**  
HEALTHCARE  
LTD

## Accessories



Permalux Hospital Pillow  
Ref. PL029



Webbed Back Rest  
Ref. 10773



Static Overbed Table  
Ref. B41 / B44 with split



Overbed Table  
Ref. A2840



Permalux Ortho Pillow  
Ref. PL030




Merlin Mattress Variator  
Ref. 10745



Hair Tray  
Ref. 5358



Medicine Chest  
Ref. 14138

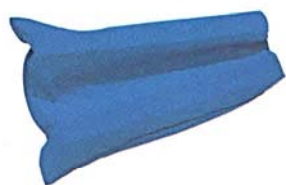

Leg Through Support  
Ref. PL032S



Lift Ease  
Ref. B54 / B52



Dressing Stick  
Ref. 10620



Tablet Splitter / Crusher  
Ref. 12912 / 14134




Permalux Protector  
Ref. PL086 Elbow / PL087 Heel



Foam Wedge  
Ref. PL089F




Button Hook  
Ref. 10630

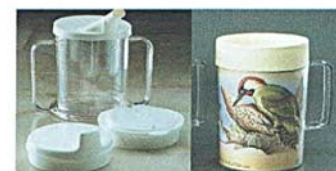

Two-Handled / Insulated Mug  
Ref. 10488 / 10458



Permalux Duvet  
Ref. PL023

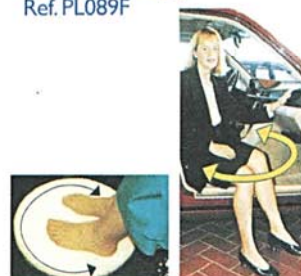

Turn Easy  
Ref. 13408



Horn / Stocking Aid  
Ref. 10604 / 10624



Extender Reacher  
Ref. 10650



Permalux Cot Bumper  
Ref. PL092BL



Slide Roll  
Ref. T004SR




Free Hand Tray  
Ref. 10490



Comfort Cutlery  
Ref. 122xx

# Bathroom

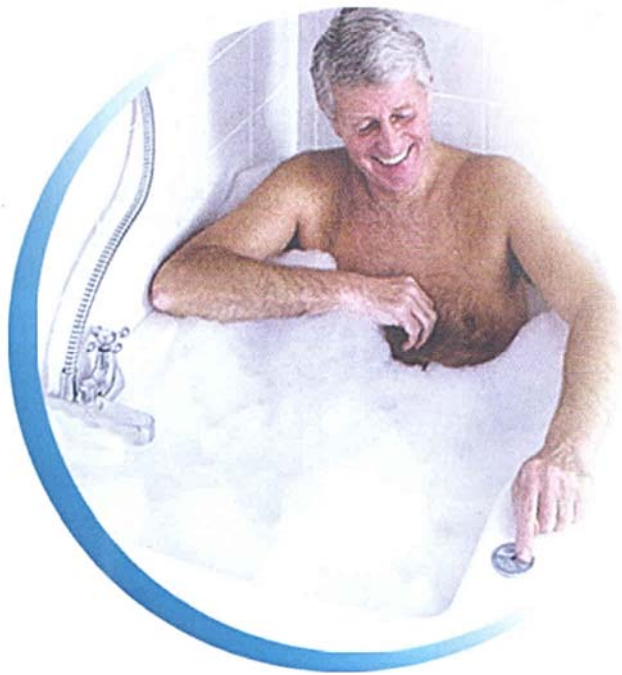

Home Mobility Lebanon

[www.homemobility.info](http://www.homemobility.info)

## Bathing




Commode



### Maxi-Drain Bath Bench

Sit, Swing legs over the bath tub walls while Sliding: Safe & Easy. Height adjustable non-slip seat with drain holes, back support, grab handle & anodized aluminium frame.

Ref. 7927A Non Padded / 7956A Padded with Commode & Pail  
O/A Dim: 71 x 47 x 83 – 93 cm  
Max. user weight: 113 Kg

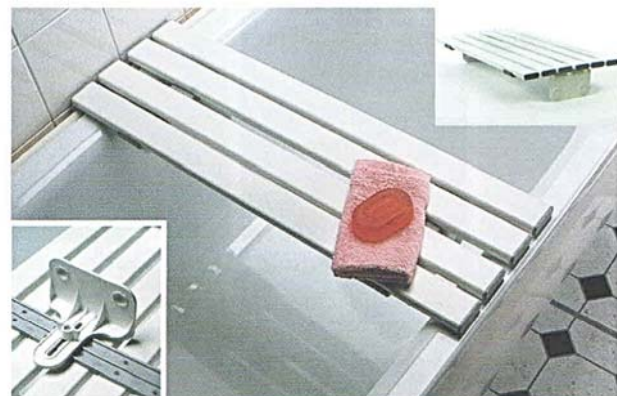

### Merlin Bathboard

Sit, Swing legs over the bath tub walls while Rotating. Made of durable reinforced plastic with slatted design for excellent drainage. Secures safely to wall tub.

Ref. 10164 / 10169  
O/A Dim: 66 x 23 x 2.5 cm / 68.6 x 38.1 x 2.5 cm  
Max. user weight: 190 Kg



### Merlin Bathseat

Secures firmly to bath floor and allows you to sit inside the bath.

Ref. 10157  
O/A Dim: 45.7 x 28 x 15.2 cm  
Max. user weight: 160 Kg



### Shower 3-Massager

Extra long reinforced tangle-free hose.

Ref. RTL12045  
Length: 213 cm



### Bath Mat

The suction cups attach securely to tub/shower floors.

Ref. 2050A  
O/A Dim: 35.5 x 91.5 cm  
Max. user weight: 150 Kg



### Turning Bath Seat

Sit, Swing legs over the bath tub walls while Swiveling. Plastic reinforced non-slip seat with excellent drainage. Strong back support. Secures safely to bath tub and locks into position.

Ref. 10174  
O/A Dim: 52 x 53 x 51 cm  
Max. user weight: 100 Kg



### Aquatec Bathboard

Sit, Swing legs over the bath tub walls while Rotating. Made of durable reinforced plastic with textured, non-slip surface. Includes blue hand grip and soap holder. Secures safely to wall tub.

Ref. 10154-10  
O/A Dim: 69.5 x 30 x 12 cm  
Max. user weight: 150 Kg




### Otter Bathing System

Bath seat with angle-adjustable back at 0°, 22.5°, 45°, 67.5° & 90°, soft fabric, height and width adjustable head cushions and removable positioning straps.

Ref. OT2000SF/OT3000SF  
O/A Dim: 43 x 107 x 65.5 cm / 43 x 127 x 81 cm  
Max. user weight: 54 Kg / 72 Kg

## Bathing

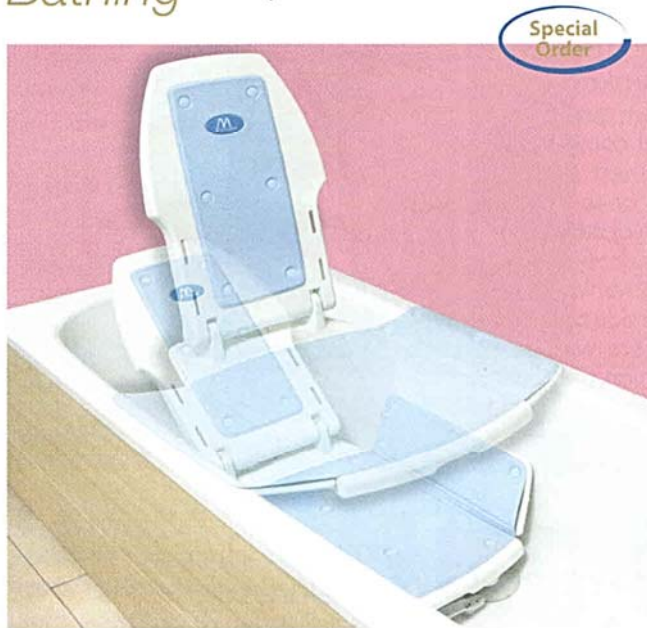

Easy to Use



Easy to Recline



Easy to Remove



Easy to Handle

### Minivator Bliss Bath Lift

If you want to enjoy the benefits of a soak in the bath without having your bathroom adapted, Bliss Bath Lift is the solution: Place it securely in the bath when needed; Raise it to upright position, leveled with the top of your bath; Sit on your comfortable anti-slip seat and move your feet and legs over into the bath; Lower the seat into the bath with the water already in it, and regain the pleasure of a relaxing bath. When you are finished, raise the bath lift to the starting position and exit the bath without having to run the water out first.




#### Ref. Bliss


Seat width: 38 cm (67.2 cm with flaps)

O/A Dim: 67.2 x 67 x 100 cm

Max. user weight: 130 Kg

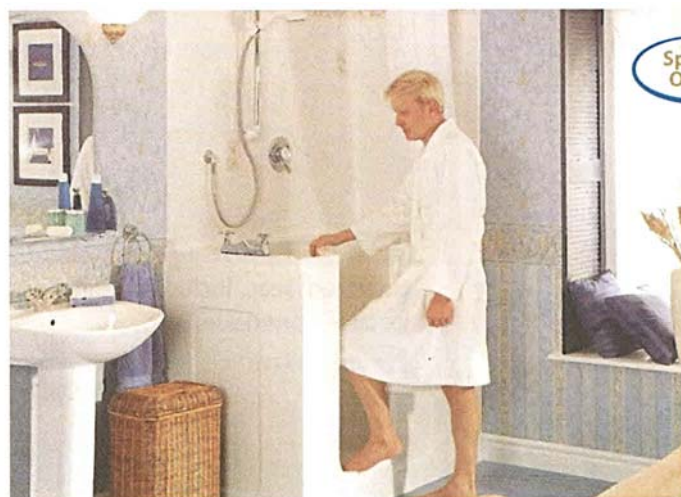

Special Order

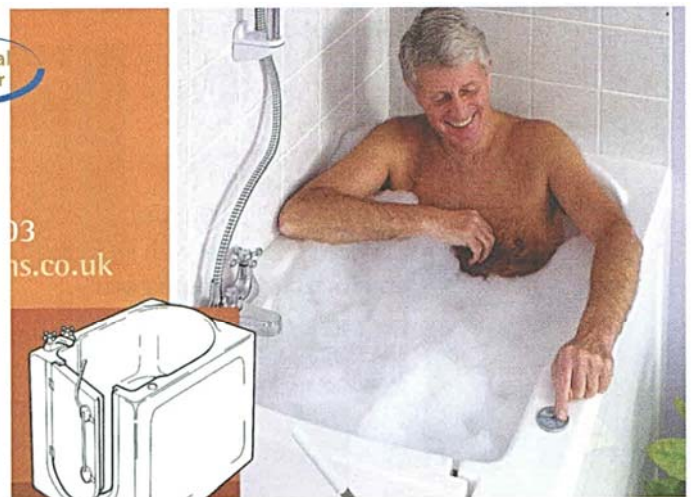

### Coniston

Walk-In Bath & enjoy a deep soak. The contoured backrest and built-in anti-slip seat allow you to remain fully supported in the upright position. Supplied with a thermostatic valve for your complete comfort and security. Optional: hydrotherapy.

#### Ref. Coniston

Bath Dim: 66 x 125 x 95 cm



جهاز حمامك  
للتجنب الوقعة



## Toileting

### Merlin Raised Toilet Frame

Stand up easily from the toilet seat: the frame is adjustable to your height and the 2 armrests give you the support you need. The plastic seat is washable and can be easily fitted.

Ref. 10380

O/A Dim: 53.5 x 42 x 60 - 80 cm  
Max. user weight: 100 Kg



### Commode Chair 4-brakes

A high quality and versatile steel commode chair; can be used over the toilet seat or as a transfer chair; comfortable moulded seat, drop-down arms and removable footrests for easy side transfer. Include 4 castors with brakes, soft cushion, pail & cover.

Ref. 10183

O/A Dim: 55 x 88 x 97 cm  
Max. user weight: 100 Kg



### 3 in 1 Steel Commode

Safety frame, commode or height adj. raised toilet seat. Contoured armrest for extra comfort. With pail & cover.

Ref. 7103A / 30216B Bariatric

O/A Dim: 55 x 44.5 x 55 - 75 cm / 65.5 x 48.5 x 64 - 81 cm  
Max. user weight: 136 / 182 Kg



### Commode Chair 2-brakes

Can be used over the toilet seat or as a transfer chair; Drop-down arms and removable footrests for easy side transfer. Include 4 castors (2 with brakes), cushion, pail & cover.

Ref. 3004/T

O/A Dim: 57 x 83 x 96 cm  
Max. user weight: 120 Kg



### Aquatec Toilet Seat

Elegant design toilet seat with 2 flip-up armrests. Made of durable, corrosion-free material.

Ref. 10098-10

O/A Dim: 60 x 43 x 31 cm  
Max. user weight: 150 Kg



### Aluminum Commode Chair

Economical lightweight commode with contoured fixed armrests. Include 4 castors (2 with brakes), pail & cover.

Ref. 2215A

O/A Dim: 51 x 49 x 86 cm  
Max. user weight: 136 Kg



### Locking Raised Toilet Seat

Economical model with wide & comfortable contoured seat

Ref. 6486R / 6487R with arms

O/A Dim: 52 x 41.5 x 14 / 29 cm  
Max. user weight: 136 Kg



### Uplift Commode Assist

Stand up or sit gently with this self-powered seat. Includes height adjustable commode frame.

Ref. CA100+7103A

O/A Dim: 55 x 44.5 x 55 - 75 cm  
Max. user weight: 100 Kg



### Male & Female Urinal

Anatomical shape; 1-liter capacity with holder.

Ref. 10844 / 10846

O/A Dim: 10 x 32 x 13.5 cm



### Commode Cushion Akros

Combination gel-foam cushion provides pressure relief while toileting.



### Bidet Bowl

Convert your toilet into a bidet. 1-liter capacity with soap dish. Simply tip back to empty.

Ref. 10838

O/A Dim: 33.5 x 35.5 x 12.5 cm

Ref. 8164561

O/A Dim: 43 x 40 x 2.54 cm  
Max. user weight: 113 Kg

## Showering



### Self-Propelled Shower/Commode Chair

4-in-1 shower / commode / transport chair or over toilet seat: Non-rust aluminum frame is height adj. with detachable flip up padded armrests & height adj. footrests. With seat & pan.


Ref. 10210C

O/A Dim: 61 x 114 x 54 - 64 cm

Max. user weight: 120 Kg

**COOPERS**  
Count on us to care

Special Order

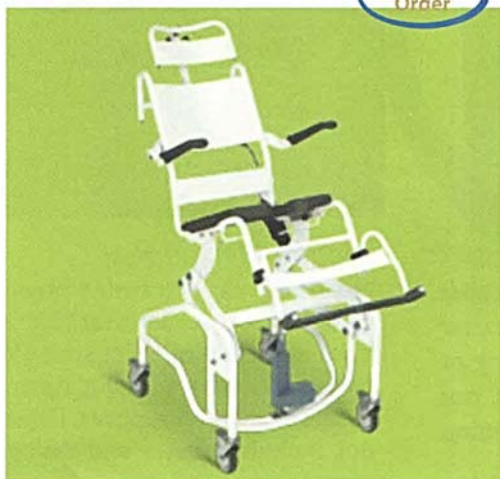

### Height Adjustable/Tilting Shower Chair

Height adjustment allows caregiver to work at a comfortable height, while 35° tilting gives correct user's positioning, maintaining eye contact during showering.

Ref. LI2161-121

Seat Dim: 45 x 45 x 58 - 104 cm

O/A Dim: 65 x 106.5 x 119.5 - 165.5 cm

Max. user weight: 150 Kg

### Shower Trolley

Adjustable back rest shower trolley, height adjustable via hydraulic pump, with HT-cover material, pvc water receiver with 1m flexible hose, 2 hinged side fences & 4 SS double braked castors; available in coated steel or stainless steel. Can be ordered as a wall-hang model, that folds when not in use.

Ref. LI2514.190

O/A Dim: 74 x 190 x 55 - 90 cm



Special Order

### Shower seat with backrest and armrests

Ergonomic glassfibre reinforced nylon shower seat with independent flip-up armrests for lateral support & easier transfer. Does not include screws and fittings.

Ref. LI2203-200

O/A Dim: 45 x 40.5 x 50 cm

Max. user weight: 150 Kg



### Padded Bath Seat

Shower comfortably with a cushioned seat and backrest, flip-up armrests; height adjustable, made of rust-resistant anodized aluminium frame.

Ref. 7945A / 7947A Commode

O/A Dim: 52.2 x 53.3 x 91 - 103 cm

Max. user weight: 136 Kg



### Shower Stool

Height adjustable

Ref. 10231C

O/A Dim: 32 cm diam.

40 - 50 cm height

Max. user weight: 113 Kg



### Turner Cubicle Seat

Wall mounted shower seat with leg support, plastic coated for corrosion protection. Folds up when not in use. Does not include screws and fitting.

Ref. 10240

O/A Dim: 46.5 x 40.7 cm

Max. user weight: 190 Kg



### Shower Chair

Light-weight height adjustable shower chair with anti-slip seat & back and built in handles. Special seat cut for personal hygiene.

Ref. 7921A / 7923A with cut

O/A Dim: 50.8 x 43 x 35.5 - 55.8 cm

Max. user weight: 158 Kg

Special Order



## Special Design



Looking for a design for your specific needs? Contact us and we'll make your life easier.

## Grab Bar

COOPERS



### Swedish Bath Rail

Use it as a safe grip to get in & out of the bath. Easy to install. Does not include screws and fitting.

Ref. 10330

O/A Dim: 37 x 18 x 60-70 cm  
Max. user weight: 100 Kg



### Toilet Frame

Surround your toilet seat with this height adjustable aluminium frame. It will help you sit or stand up.

Ref. 10397

O/A Dim: 53.5 x 40 x 61 - 70 cm  
Max. user weight: 160 Kg



### Wall to Floor Grab

Fit it next to the toilet seat. It will give the user support when sitting or standing up. Does not include screws and fitting.

Ref. 10180

O/A Dim: 8 x 75.5 x 74 cm  
Max. user weight: 100 Kg




### Looped Hinged Rail

Need an easy transfer from wheelchair to toilet seat? Fit 2 from each side of toilet seat, lift one up, transfer & pull it down to use as a lateral support. Does not include screws and fitting.

Ref. 10328

O/A Dim: 10.2 x 74 x 22.7 cm

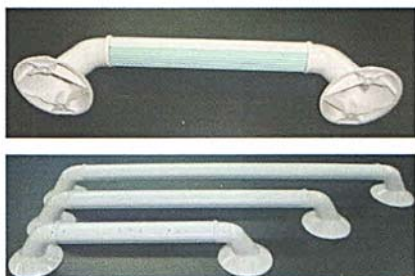

### Cefndy De-luxe Plastic Grab Bar

Fit it anywhere in your toilet, at entrance, in bath, near lavatory. Provided with a "soft-feel" ribbed underside grip. Does not include screws and fitting.

Ref. G26 / G27 / G28

O/A Dim: 4 x 30/45/61 x 8 cm  
Max. user weight: 157/125/100 Kg



### Cefndy Modular System

Secure your path walk. Add-on flexible grab bar system that can fit to any pathway length. Does not include screws and fitting.

Ref. 72xxx

O/A Dim: 4 x 8 cm  
Max. user weight: 100 - 157 Kg



### L-Shaped grab bar

Fix it on your shower wall. It will provide you a big surface to hold on. Made of heavy duty steel construction with rust resistant enamel coating. Does not include screws and fitting.

Ref. 4063A R / 4073A L

O/A Dim: 6.5 x 40 x 82 cm  
Max. user weight: 113 Kg

Access 



Home Mobility Lebanon



Stairlift Lebanon by Handicare

[www.homemobility.info](http://www.homemobility.info)

## Hoist



### Sunlift Patient Lift

Transfer your beloved in a safe and dignified way without back injury, from bed to chair & vice versa. You can even pick him up from the floor. Battery operated with emergency stop button.

Ref. G130E/G150E/G175E

O/A Dim: 57.5-105 x 112 x 125.5-186.5 cm

66.5-120.5 x 121.5 x 125.5-193 cm

66.5-122 x 134 x 142.5-198.5 cm

Max. user weight: 130/150/175 Kg

**Joerns Healthcare**



### Oxford Slings



### Voyager Ceiling Hoist

Can be portable & mounts easily on posts, or can be fixed on a ceiling track. The wide range of track options covers all transfer needs.

Ref. Voy-9800/420/550/800



Max. user weight: 200/190/250/360 Kg



### Easy Track



### Ceiling Track

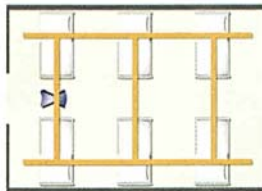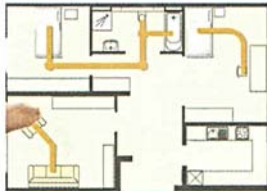

### Emergency Stop



### Safety Spreader Bar



### Push Handle



### Legs Adjustment



### Floor Pickup

### SAFE & EASY TRANSFER



صار فيك تنقل الشخص من  
السريير للكرسي بطريقة سهلة  
وأمنة وبدون ما يبرق ظهره!



**SUNRISE MEDICAL**

## Hoist

**Veziris**  
HEALTHCARE

**Joerns**  
Healthcare



Special Order

### Body Up


New lift concept, transfer from/to bed, chair, toilet, shower and car.

Ref. B1346

O/A Dim: 61 - 73 x 88 x 99 cm

Max. user weight: 100 Kg

New

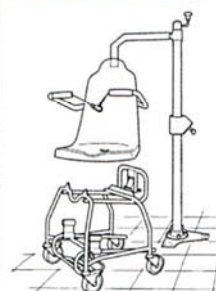

### Oxford Dipper Pool Hoist

Get you out of the pool in a safe and dignified way. Integrate fully with pool environment with 3 floor fixing options that can fit in almost all locations. Use it with a sling or Ranger transporter.



Pool

Special Order



Sea

### JoB

Now you can swim without worrying how to get in & out of the pool or the sea: J.O.B. is the original Beach Wheelchair, made to move disabled and elderly people in every situation: from sand to water, from snow to off-road. JoB easily folds up so it can be transported in small vehicles.

Ref. Job

O/A Dim: 93 x 107 x 93 cm

Max. user weight: 100 Kg

**NEATECH**  
SOLUZIONI PER LA MOBILITÀ



New



Special Order

### I-Swim

The battery powered lift will get you in & out of the pool in an easy and safe way.


Ref. i-swim

O/A Dim: 64 x 125 x 151 cm

Max. user weight: 140 Kg

**AUTOLIFT**

## Ramp

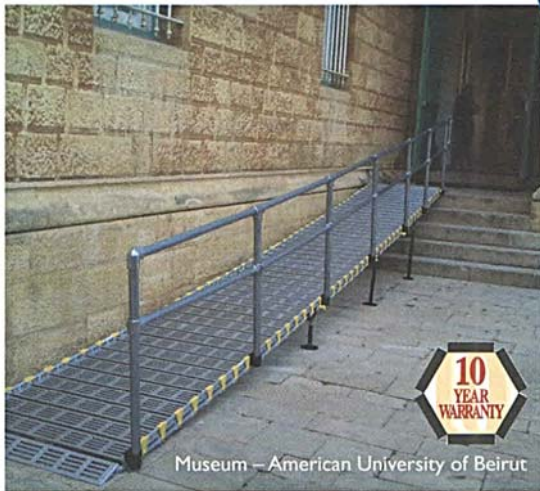


Museum – American University of Beirut

**10  
YEAR  
WARRANTY**

Special  
Order



Special  
Order

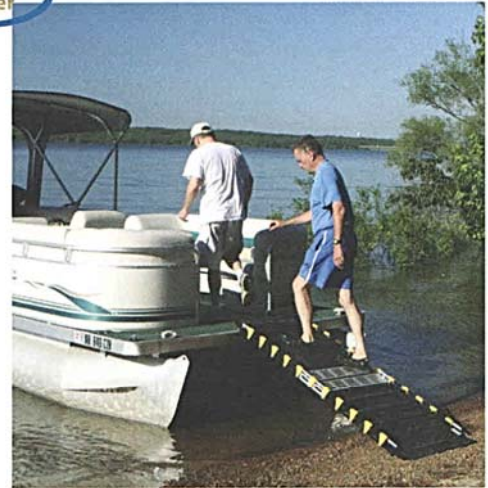

### Roll a Ramp

Customize this unique ramp for each individual circumstance by simply adding on it a link of the length you want. After usage, you can split it, fold it & transport it in your car. Constructed of aircraft quality anodized aluminium, it is extremely durable and has 10 years warranty.

Ref. 12" pair/26"/30"

Max. user weight: 454 to 908 Kg

Boat



Special  
Order



New



Piuma the feather light all aluminium lift



Special  
Order

### MiniFold

Install it in your existing van and spare the cost of a lowered floor conversion. All you need is a min. door clearance of 72.4 cm and a fairly flat floor surface inside door. Ramp can unfold automatically or manually.

Ref. Power / Manual MiniFold

Max. user weight: 454 Kg



### Piuma

The feather light all aluminium lift

Ref. BA11401

Platform size: 90 x 140 cm

O/A Dim: 127 x 142 x 152 cm

Max. user weight: 350 Kg

ما تخلي العجز يوقف حياتك.  
التغيير بمحيطك بسهل عيشتك!



## Straight Stair Lift



### Simplicity 950 Indoor

Enjoy the very essence of a stairlift without compromising on quality or style.



Ref. M950

Seat: 48 x 42 x 39-49 / Swivel: 68 / Folded: 28.5 cm

O/A Dim: 62.5 x 58 x 90.5-105.5 cm (WxDxH)

Speed: 0.12 m/sec. Stair Angle <56°

Max. user weight: 140 Kg



Special Order

### Minivator 1000 Indoor / Outdoor

Now you can install your stairlift outside! The fully waterproofed unit comes with 1 pendant and cover to help keep it clean.

Ref. M1000

Seat: 44 x 40 x 45.5-53 / Swivel: 72.5 / Folded: 40 cm

O/A Dim: 56 x 68 x 162.5-170.5 cm (WxDxH)

Speed: 0.12 m/sec. Stair Angle <56°

Max. user weight: 137 Kg. (Opt 200 Kg)



Your stairs are no longer a problem! Just sit on the seat, put your seat belt, press the toggle & enjoy the ride. When you reach the top of the stairs, push the lever available on both sides, swivel & dismount safely. When not in use, you can fold it up. You can also operate it via remote control. All Handicare stairlifts have a battery back-up in case of power failure.



Unique Toggle



Swivel Seat



Safety Belt



Safety Edge




Key Switch



Pendant



Hinge Track

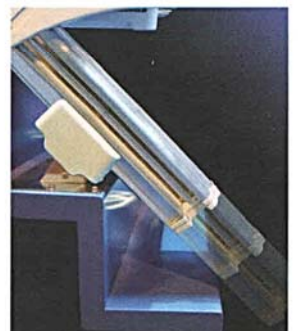

Slide Track

صار في كرسي  
تطلعك الدرج!



## Curved Stair Lift



Designed with safety as a standard: the seat is upholstered with flame retardant wipe clean fabric; the footplate is covered with anti-slip hard wearing carpet and is equipped with safety edges that will stop the lift should it meet an obstacle; the hide-away key switch prevents unauthorized use; the self-diagnostic intelligent system will easily identify issue leading to a quicker & more efficient rectification.



Manual Swivel




Unique Toggle




Powered Swivel Seat



Simple to Store

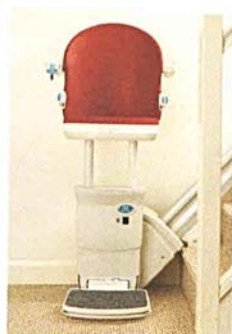

Perch Seat

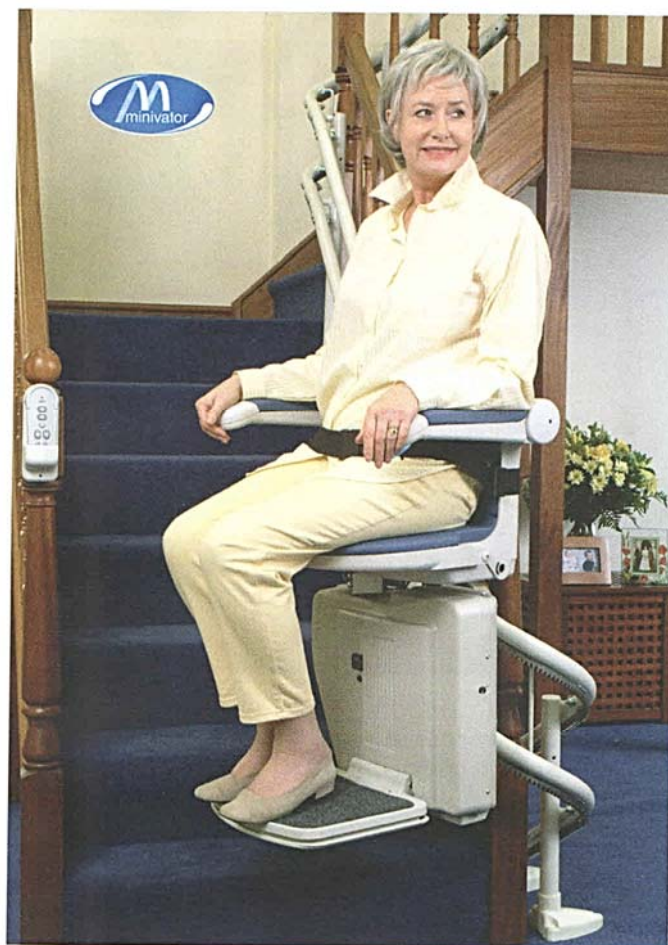

Minivator 2000 Indoor / Outdoor

Best of its class stairlift that can be fitted in locations that no other manufacturer can accommodate, thank's to Minivator skilled & dedicated research team.

### Seat Color Options



Cocoa



Sapphire



Ruby



Slate



Sand



Jade

### Ref. M2000


Seat: 44 x 40 x 50.9-58.4 / Swivel: 73.7 / Folded: 40 cm

O/A Dim: 56 x 66.8 x 100-108 cm (WxDxH)

Speed: 0.12 m/sec. Stair Angle <56°

Max. user weight: 115 Kg. (Opt 137 Kg)

### Rail Options

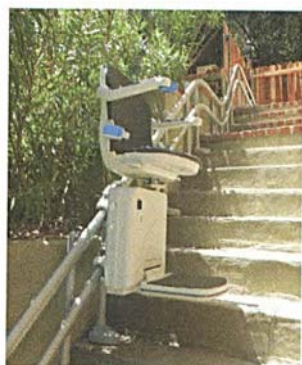

Outdoor



90° Configuration



Spiral Configuration



Choose Your Rail Color



180° Configuration



Top Overrun



Hinged Rail

## Curved Stair Lift



**New** Up to total rail length of 35 meters



Freelift's stairlifts have an attractive single tubular rail and a wide variety of seats, in a choice of materials and colors. A child's seat is also available.

### Ref. Freelift

Minimum Stair Width: 69 cm

Max. user weight: 125 Kg



45° Downward Facing Seat

### Seat Options



Classic



Comfort



Elegance



Panda


### Rail Options



Hinged Rail



90° top parking bend

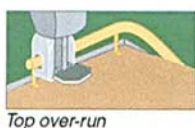

Top over-run



Standard



Drop Nose



Bottom over-run



90° parking bend bottom



180° parking bend bottom



Rail Color



Made in Britain

Special Order

**Joystick**  
This easy to use control can be positioned on left or right arm.

**Key Switch**  
The key can be removed to prevent unauthorized use.



**Space Saving Design**  
The armrests & foot support can be folded away easily, clearing space on the stairs.



**Safety Belt**  
Provided for your protection & security.



**Seat Belt Extension**  
Optional.



**Manual Swivel**  
Push down lever allows seat to be swiveled enabling safe mounting & dismounting.

**Powered Swivel**  
Optional, operated by joystick.

**Footrest Linkage**  
Optional, folds up with seat.

**Off-centre Footrest**  
Optional unique design for tight steps.

## Inclined Platform Lift



**Ref.V64**

Platform: 83x70 / 105x77 / 125x80 (WxD)

Folded: 37 cm

Speed: 5 m/min. Stair Angle <50°

Max. user weight: 200 - 230 Kg



Special  
Order



### V64 / V65

Inclined platform lift is an ideal solution for wheelchair users to enjoy an adequate social life by allowing them to get up and down stairs on their own. Simply roll in the wheelchair on the platform, pull down the arm & push the button. Upon arrival at stair end, pull up the arm and roll the wheelchair out. It can be installed in already existing structure or outdoor.

**Ref.V65**

Platform: 83x70 / 105x77 / 125x80 (WxD)

Folded: 37 cm

Speed: 8 m/min. Stair Angle <50°

Max. user weight: 200 - 230 Kg

## Stair Climber

Special  
Order



## Vertical Platform Lift

### Silver

Silver allows trouble-free access, with its stainless steel structure design and state-of-the-art active and passive safety features thanks to the well proven Vimec mechanical drive systems. Available in 80 cm and 100 cm lift range and does not require a pit for installation.

**Ref. Silver**

Platform Dim: min: 120 x 85 cm

Max: 140 x 125 cm

Speed: 0.04 m/sec.

Max. user weight: 250 Kg



### Roby

Hook the wheelchair on Roby and climb the stairs effortlessly. Innovative design makes it extremely comfortable, due to its constant speed, safe and easy to operate. You can easily dismantle it and take it in your car.

**Ref.T09**

Stair gradient < 35°, width > 82, landing depth > 97 cm

Weight: 26 Kg drive unit, 9 Kg batteries, 10 Kg steering bar

Speed: constant 5m/min.

Max. capacity: 130 Kg

### Steppy

Steppy is the ideal way of assuring access to homes or businesses with lifting ranges of up to 2 meters. Full platform stability and a smooth ride are ensured by the cushioned arrival. The user's safety is guaranteed by all the latest active and passive safety systems.



Special  
Order



**Ref. Steppy**

Platform Dim: min: 126 x 81 cm

Max: 220 x 91 cm


Speed: 0.06 m/sec.

Pit: min: 22, max: 36 cm

Max. user weight: 400 Kg



## Elevator

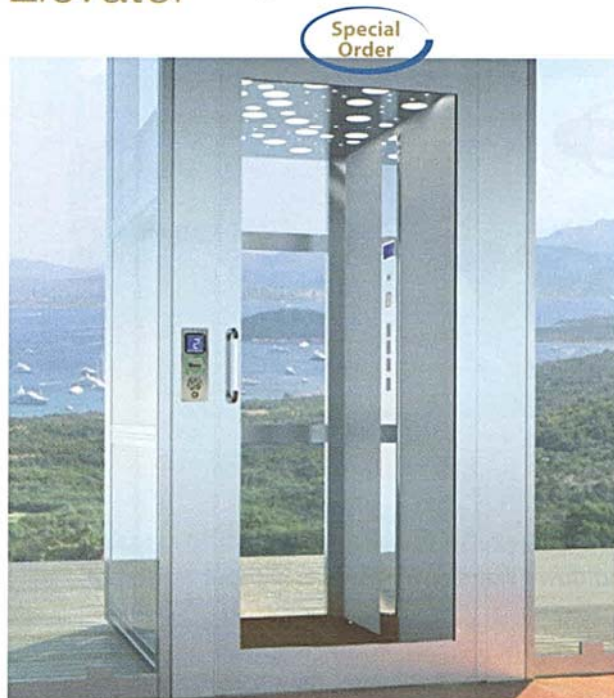

## Cabin Doors




## Cabin Platforms



## Self Supporting Structure



## Cabin Walls

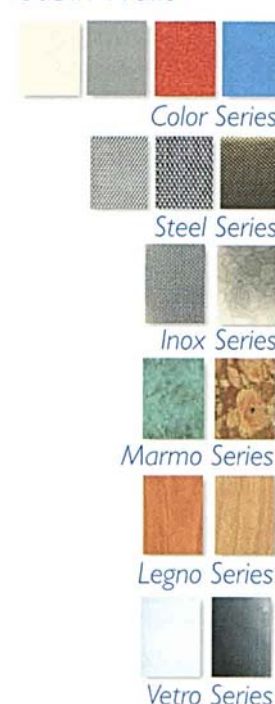

## Easy Move

Easy Move hydraulic lift fits in tight spaces with simple tailor-made installations, both indoor and outdoor. The Self Supporting Structure resolves the need for special elevator room or masonry works, can go up to 5 floors with a pit of just 14 cm: it is an ideal solution for renovated & new buildings.

Customize cabin size & colors to match your décor with the world-famous Italian design. Quiet and discreet, it can be used in complete comfort and safety thanks to the state-of-the-art technical features that safeguard users in all stages of use, with protected, user-friendly controls. More than one thousand installations made worldwide.

## Ref. Easy Move

Platform dimensions: min: 90 x 90 cm, max: 120 x 140 cm

Maximum travel: 5 stops (14.6 m)

Speed: 0.15 m/sec

Max. user weight: 400 Kg



**VIMEC**  
Easy moving



Control Box

## EcoVimec

EcoVimec, your ecological lift, operates without hydraulic oil, using manufacturing process with waste material for recycling & low-energy consumption technology.

## Ref. EcoVimec

Platform dimensions: min: 90 x 90 cm, max: 110 x 140 cm

Maximum travel: 5 stops (14.6 m)

Speed: 0.15 m/sec

Max. user weight: 400 Kg

## Regain Your Freedom

**New**



**TopChair**



**Special Order**



### TOPCHAIR-S

The only powered wheelchair in the World capable of going up and down stairs and sidewalks without assistance.

Ref. TSC-52

O/A Dim: 70 x 115 x 107 cm

Max. user weight: 110 Kg

**New**

### QuadriX

Seeking for thrills ? Need to escape ? Or simply enjoying the pleasure of discovering natural spaces ? Either on your own, or with family or friends, QuadriX will give you the opportunity to enjoy the thrills of 4-wheel off-road downhill mountain biking..

Ref. Watts

**QUADRI X**

**Special Order**



**Special Order**

**New**



**Ferriol Matrat**

### Joelette

is a single-wheel off-road wheelchair that allows the practice of hiking for all disabled, even heavily dependent, with the help of two companions.

Ref. C10P01

O/A Dim: cm

Max. user weight: Kg

*Mobility* 



Home Mobility Lebanon



Quickie Wheelchair Lebanon

[www.homemobility.info](http://www.homemobility.info)

## Walking Aid



Slow Down Brakes



Brake Extension

Rental Available



One Hand Brake



Ref. 12070  
Seat height: 63 cm  
O/A Dim: 59 x 74.5 x 78-93.5 cm  
Max. user weight: 125 Kg



### Dolomite Wheeled Walker

Get the best upright posture & a reliable walking support with this folding walker, recognised world-wide for its ergonomic design and safe breaking system. Transport your things in the basket. Met some friends on the way? Pull the brakes, fold the seat & sit: your back is still supported.



**Shoppy Wheeled Walker**  
Economical version of Dolomite wheeled walker without back support.

Ref. 12300  
Seat height: 58 cm  
O/A Dim: 59 x 66 x 81-99 cm  
Max. user weight: 126 Kg



603630A wheels

Rental Available

**Folding Walker**  
Durable, lightweight aluminum H-frame, safe hand hold, dual-release, height adjustable.

Ref. 616070A Adult / 616071A Jr  
O/A Dim: 54/48 x 40/34 x 84-101/69-88 cm  
User height: 160-188 / 140-160  
Max. user weight: 136 Kg



**Imperial Walker**  
Extra-wide, dual-release, height adjustable.

Ref. 604070A  
O/A Dim: 58.5 x 45 x 81-99 cm  
Max. user weight: 272 Kg



**Kid's Walker**  
Folding, height adjustable with 2 front wheels.

Ref. A135  
O/A Dim: 50 x 57 x 54-78 cm  
Max. user weight: Kg



**Fisher Support Stick**  
Specially designed handle that spreads the pressure evenly across the palm of the hand. Height adjustable.

Ref. 7655 L/7657 R  
O/A Dim: 7 x 13 x 73-94 cm  
Max. user weight: 127 Kg



**Foam Swan Cane**  
Coopers walking stick with foam handle, height adjustable.

Ref. 7674  
O/A Dim: 3.5 x 16 x 74.5-100 cm  
Max. user weight: 100 Kg



**Wood Handle Cane**  
Coopers walking stick with wooden handle, height adjustable.

Ref. 7665  
O/A Dim: 2.2 x 13 x 72-98 cm  
Max. user weight: 100 Kg



**Tetrapode**  
Self-standing, provide a stable support, right, left & height adjustable.

Ref. 7808  
O/A Dim: 17 x 25 x 72-97 cm  
Max. user weight: 127 Kg

COOPERS

Ferrule  
Ref. 1.5



**Aluminum Cane**  
Contoured vinyl handle provide secure hand hold. Height adjustable.

Ref. 6220A  
O/A Dim: 3 x 14 x 76.5-98 cm  
Max. user weight: 113 Kg



**Universal Aluminum Crutches**

Ref. 3610M / 3612 youth  
Max. user weight: 113Kg



**Coopers Comfy**

Ref. 8351  
O/A Dim: 10x18x97-119 cm  
Max. user weight: 127 Kg

**Elbow Crutches**

Ref. 128 / 130  
O/A Dim: 13x20x87-116 cm  
Max. user weight: 100 Kg



## Standard Wheelchair

**Breezy**  
Premium products to Quality service

Breezy 250/300 folding wheelchairs have removable flip-up armrests, swing-in/out footrests.



**Breezy 300 Aluminum Chair**  
Simply light with a modern look.

Ref. B305 / B315 half-fold back  
Seat: 37/40/43/46/49/52 x 41 cm  
O/A Dim: 34 cm (folded) x 106 x 95 cm, 13 Kg  
Max. user weight: 120 Kg



**Breezy 250 Steel Chair**  
New design, strong & durable.

Ref. B250  
Seat: 37/40/43/46/49/52 x 42 cm  
O/A Dim: 34 cm (folded) x 103 x 94 cm, 16.5 Kg  
Max. user weight: 125 Kg



**Breezy 250R Steel Reclining Back**  
Includes elevating legrests.

Ref. B250R  
Seat: 37/40/43/46/49/52 x 42 cm  
O/A Dim: 39 cm (folded) x 103 x 128 cm, 18.5 Kg  
Max. user weight: 125 Kg



**Breezy 90 Steel Chair**  
Removable flip-up arms & swing-out legrest.

Ref. B90  
Seat: 40/42/45 x 40 cm  
O/A Dim: 28 cm (folded) x 97 x 91.5 cm, 15.8 Kg  
Max. user weight: 125 Kg



**Companion Chair**  
Folding back, removable SA footrest.

Ref. EJ766-1  
Seat: 46 x 40 cm  
O/A Dim: 51 cm x 91 x 100 cm, 11 Kg  
Max. user weight: 113 Kg

## Positioning Wheelchair

Special Order



**Breezy Relax2**  
Comfortable with adjustable seat height/depth, back height & length arm-pads.

Ref. Relax2  
Seat/Weight: 41-53 x 42-50 cm / 34 Kg  
Tilt/Recline: -3° -22° / 0°-25°  
Max. user weight: 145/170 Kg



**SR45**  
Intelligent Rotation in Space technology means the seat frame rotates around the user's center of gravity.

Ref. EIZ9  
Seat/Weight: 36-53 x 36-56 cm / 17 Kg  
Tilt/Recline: 45° / 5°-30°  
Max. user weight: 120 Kg




**Quickie IRIS**  
Portable & lightest Tilt-in-Space chair with 40° or 55° I.R.I.S..

Ref. EIZ4-2  
Seat/Weight: 36-56 x 36-56 cm / 17 Kg  
Tilt/Recline: 40°/55° / 5°-30°  
Max. user weight: 113 / 158 Kg

## Sport Wheelchair

Special Order

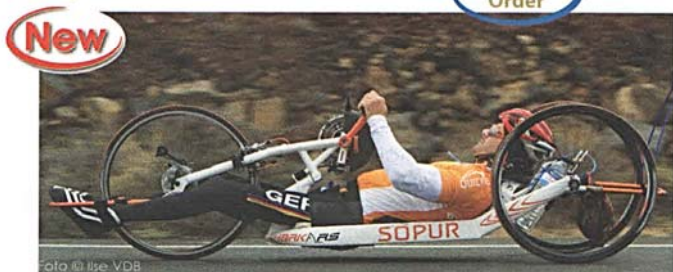

**Shark RS**



**Match Point**



**All-Court**

**QUICKIE**  
Live without limits.

**SUNRISE MEDICAL**

## Active Rigid Wheelchair

Special Order

QUICKIE  
Live without limits.

As Individual  
as Your  
Fingerprint

SOPUR.  
Live without limits.

CRASH  
TESTED



reddot award  
product design



### Helium Gen 2013

Minimalist design and impressive driving performance: extremely strong and sturdy.

Ref. Helium  
Seat/Weight: 32-46 x 34-48 cm / > 6.4 Kg  
Max. user weight: 125 Kg



### Quickie Q7 NextGen

Designed for superior handling, efficiency, and speed, Q7 is the lightest wheelchair we've ever manufactured.

Ref. EIR4  
Seat/Weight: 30-50 x 30-50 cm / > 6 Kg  
Max. user weight: 120 Kg



### Argon 2

The ideal combination of a rigid performing and multi-adjustable wheelchair: the wheelchair that's adaptable to you!

Ref. Argon  
Seat/Weight: 30-50 x 30-50 cm / > 9 Kg  
Max. user weight: 125/140 Kg




### GP/IGP

Cruise in a classic Quickie, experience its durability & comfortable ride quality.

Ref. EIGP  
Seat/Weight: 35-50 x 35-50 cm / > 11 Kg  
Max. user weight: 113 Kg



Improper positioning in a wheelchair can cause any of the postures described.

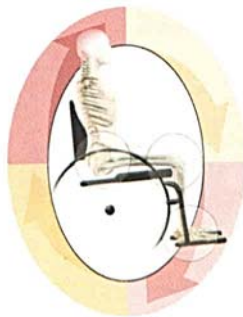

### Life R

Modular design, adjustable seat depth makes it perfect for growing teenagers.

Ref. LifeR  
Seat/Weight: 34-48 x 36-50 cm / > 10.7 Kg  
Max. user weight: 125 Kg

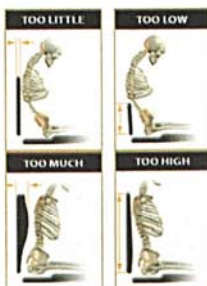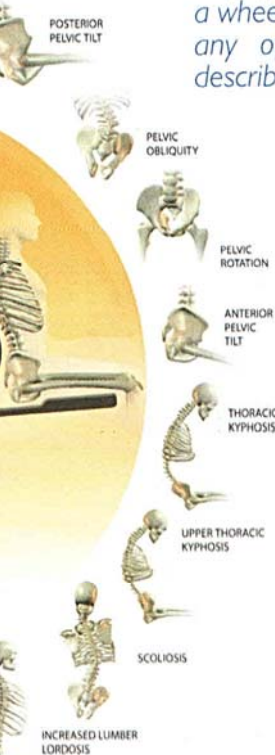

خلينا نساعدك لتختار الكرسي المناسب لحالتك. ومن صنعك بها على قياسك!



## Active Folding Wheelchair



### Xenon

The lightest folding wheelchair in the world.

Ref. Xenon

Seat/Weight: 32-46 x 34-50 cm / > 8.4 Kg  
Max. user weight: 125 Kg



### Quickie 2

A tried and tested classic, available in an endless number of configurations.

Ref. EIQ2

Seat/Weight: 30-56 x 30-50 cm / > 11.5 Kg  
Max. user weight: 120 Kg



### Neon

Small. Smaller. Smallest! Able to fit practically anywhere.

Ref. Neon

Seat/Weight: 30-50 x 30-50 cm / > 9.7 Kg  
Max. user weight: 115/140 Kg



### Revolution

The rigid chair that folds! Goes even inside airplanes.

Ref. EI2D

Seat/Weight: 35-50 x 35-45 cm / > 12.5 Kg  
Max. user weight: 113 Kg



**New**

### Life

Designed to meet the needs of modern life. Can be configured to suit you & your life style.

Ref. Life

Seat/Weight: 36-50 x 36-50 cm / > 12.1 Kg  
Max. user weight: 125 Kg



**New**

### Easy 200

The Benchmark in active mobility ... made even more user-friendly.

Ref. Easy200

Seat/Weight: 32-54 x 34-50 cm / > 12.3 Kg  
Max. user weight: 140 Kg



**New**

### HeliX 2

The configurable active chair designed for efficient fleet use.!

Ref. Helix 2

Seat/Weight: 38-52 x 41-51 cm / > 15.1 Kg  
Max. user weight: 125 Kg



### Quickie M6

Quickie's customization with heavy duty strength for bariatric users.

Ref. EIFI

Seat/Weight: 56-76 x 45-56 cm / > 24 Kg  
Max. user weight: 295 Kg

## Kid's Wheelchair

Special Order

zippie.  
Explore your world!

SUNRISE  
MEDICAL

CRASH TESTED



### Salsa M

Power mid wheel chair with outstanding indoor maneuverability.

Ref. Salsa M

Seat/Weight: 30-40 x 30-40 cm / 113 Kg  
Tilt/Speed: -5° - 25° or 0° - 30° / 6 & 10 kph  
Max. user weight: 75 Kg

### TS

Zippie's Original pediatric tilt-in-space wheelchair.

Ref. EIZIA

Seat/Weight: 25-45 x 33-50 cm / > 13.1 Kg  
Max. user weight: 75 Kg

### Levia Basculante Kid

The dynamic chair that absorbs spastic movement, grows with your kid & comfortably supports his positioning needs.

Ref. S081LBKE

Seat/Weight: 30-38 x 30-42 cm / > 15 Kg  
Max. user weight: 80 Kg



### Zone

Get in the Zone with the lightest Zippie designed for active kids and teens.

Ref. EIRZI

Seat/Weight: 20-63 x 20-40 cm / > 6.3 Kg  
Max. user weight: 75 Kg

### Kidz

Explore with the smallest Zippie perfect for independent toddlers.

Ref. EIKZ

Seat/Weight: 25-35 x 20-35 cm / > 6.3 Kg  
Max. user weight: 68 Kg

### Youngster 3

The lightest folding wheelchair with growth adjustability.

Ref. Youngster

Seat/Weight: 22-40 x 24-40 cm / > 8.5 Kg  
Max. user weight: 85 Kg



### Jay Fit System / Zip cushion

Configure the best seating based on your child's needs.

### Kid Kart Xpress

Early Intervention mobility for toddlers up to 25 Kg.

Ref. EIXP

Seat/Weight: 15-31 x 14-28 cm / > 15.4 Kg  
Max. user weight: 25 Kg

### Mighty Lite

Light weight and portable mobility device.

Ref. EIZ15

Seat/Weight: 30-44 x 29-43 cm / > 13.1 Kg  
Max. user weight: 68 Kg

## Power Wheelchair

Special Order

QUICKIE  
Live without limits.

SUNRISE  
MEDICAL

New



### Jive M

The ultimate driving experience with SpiderTrac suspension and power tilt, recline, ELR & seat height. Also available in Rear & Front wheel and FXL versions, up to 240 Kg.

Ref. JiveM

Seat: 40-56 x 40-56 cm / Height 42-48 cm  
Weight > 115 Kg  
Speed: 6, 10, 13 kph / Radius 56 cm  
Max. user weight: 160 Kg



New



### Salsa R2

Compact for indoors, powerful for outdoors.

Ref. SalsaR2

Seat: 41-51 x 41-51 cm  
Speed: 6, 10 kph  
Radius 120 cm  
Max. user weight: 140 Kg



### P-222SE

Rule your world with speed and style: rear-wheel drive, sport frame and manual conversion kit.

Ref. EIPW2

Seat: 35-61 x 35-51 cm  
Speed: 13 kph / Radius 63.5 cm  
Max. user weight: 158 Kg



### S6 Series

Rear-wheel drive powerful performance and superior seating with an edge: power tilt, recline & ELR.

Ref. EIS62

Seat: 35-61 x 35-56 cm  
Weight: 58 Kg  
Speed: 8, 10, 13 kph  
Radius 61 cm  
Max. user weight: 136 Kg



Rental Available

### F35

Enjoy life indoor & outdoor in an economical way. Width & height can be adj. on site.

Ref. F35

O/A Dim: 60 x 106 x 95 cm - 63 Kg  
Speed: 6 Kph / Radius 90 cm  
Max. user weight: 120 Kg



New

Compass Sport  
Center-wheel drive with captain's high-back seat.

Ref. GP605CC

Seat: 40-50 x 40-50 cm  
Weight > 41.7 Kg  
Speed: 6.4 kph  
Radius 49.5 cm  
Max. user weight: 136 Kg



### LiteRider ptc

Powered Personal Transport Chair with stadium style seat, quickly disassembles.

Ref. GP162

Seat: 43 x 40 cm  
Weight > 30 Kg  
Speed: 5.6 kph  
Radius 66 cm  
Max. user weight: 136 Kg

New



### Buzzaround Lite scooter

Features our easy "Pop! Pull! Push!" wireless disassembly for transport.

Ref. GB106

Seat: 43 x 40 cm  
Weight > 44 Kg  
Speed: 6.4 kph  
Radius 79 cm  
Max. user weight: 113 Kg



GOLDEN  
TECHNOLOGIES

## Cushions & Backs



Special Order

### J3 cushion

Anatomically designed with the Science of Seating, J3 features outstanding pressure distribution, reduced surface tension, heat & moisture dissipation & optimal stability.



### Jay® Easy™

Moderate positioning & protection.



New

### Jay Union

Jay flow fluid and memory foam for high performance skin protection and positioning.

Pre-contouring


Positioning

Capability  
Tension Reduction

Microclimatic  
Material

Segmentation

Hydrostatic Force

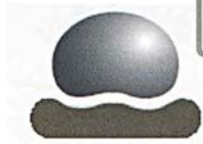

Provides pressure redistribution, positioning and stability



Provides proper positioning for postural needs



Tension-reducing materials maintain superior pressure reduction



Allows the passage of air, heat & water from the user



Help localize change to minimize instability and the spread of tension



Equalizes loading forces to reduce peak pressures.




### Jay® J3™ Back

Special Order

Spine-Align your back: set the proper angle, shape & fill in the gaps with foam. The revolutionary Jay mount hardware adjust to all posts thanks to its innovative clamp.

الأعدة بالكرسي مهم ان تكون جالسة وتأمّنك الحماية اللازمة من العقر!

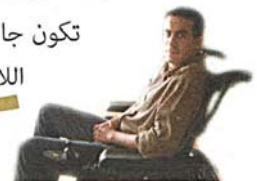

## Easy Stand



Special Order

### Glider


State-of-the-art active standing technology, designed for easier transfers & adjustability.

Ref. PNG50187

Size: Medium, Large

Max. user weight: 90, 127 Kg

### The Medical Benefit of Standing

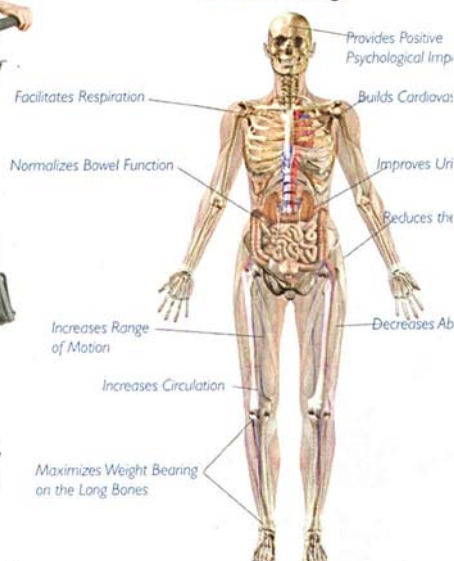

EasyStand

New



Special Order

### Bantam

The only standing frame to combine a sit-to-stand transition with the added function of an supine stander.

Ref. PY5500

Size: Extra Small, Small & Medium

Max. user weight: 23, 45, 90 Kg



*Check Our Latest News*  
**[www.homemobility.info](http://www.homemobility.info)**



**Home Mobility Lebanon**

**home  
mobility®**  
*Making Life easier*

# Rental



## Home Mobility Rental Program:

- ✓ Wide Range of High End Quality Products
- ✓ Checked, Cleaned & Disinfected after Each Use
- ✓ Meet the Highest Safety Standards
- ✓ Personalized Solutions for Your Specific Needs



Rental at Home Mobility  
*The best care  
at home...*