

The Seats

Classic

The best-selling Classic seat has an attractive leather grain effect finish. Made from durable, non-slip polyurethane, it is hard wearing and easy to keep clean.



Elegance

The new Elegance seat is lightly padded and offers a harmonious combination of elegant design and sophisticated comfort.



Comfort

The ergonomically designed Comfort Chair™ has a self lifting feature that enables the user to be easily lifted to their feet, or safely lowered into a seated position.




Panda

The Panda seat is available in four sizes, suitable for children from 1 to 9 years old.



N.B. The weight limit for the above seats is 125kg (19.5 st).

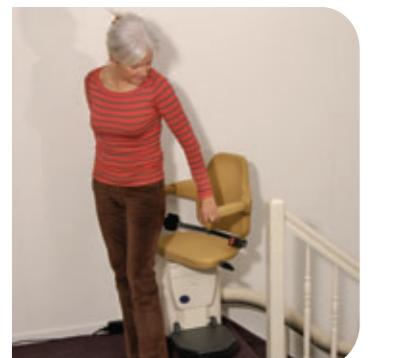

With a single movement you turn the seat and stand up easily and safely, at exactly the right place.



The footrest and carriage are fitted with safety edges so that the stairlift will automatically stop if it comes into contact with any obstacle on the stairs.

Also Available From Handicare

MOVING AND HANDLING



BATHING SOLUTIONS



STRAIGHT STAIRLIFTS



MOBILITY



Freelift Stairlifts



82, First Avenue
Pensnett Estate
Kingswinford
West Midlands
DY6 7FJ




Textphone (180015) 020 7490 8913 ext 05
Voicephone (18002) 020 7490 8913 ext 05
Dealer Support 0844 225 3121
Enquiries@handicare.co.uk
User Support 08700 118282
www.handicare.co.uk




Please be aware that actual colours may vary slightly from those shown. All information is provided in good faith but does not constitute a guarantee. Handicare reserves the right to alter designs and specifications without prior notice.

The Rails

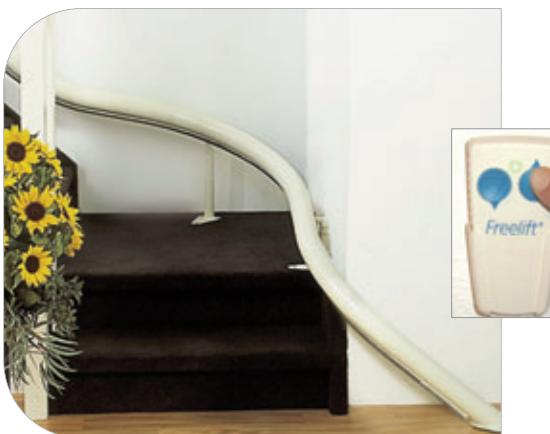

Van Gogh

The Van Gogh stairlift is designed for use on the outside bend of curved staircases. It consists of a single tubular rail, which is less than 8cm in diameter, and provides a safe, smooth and comfortable ride up even the most complicated of staircases. It runs almost silently and the compact design allows plenty of space for other stair users.

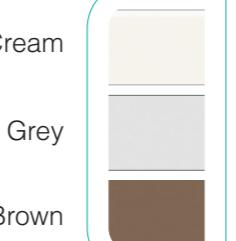


Rembrandt

The Rembrandt stairlift has the same luxurious and low-noise qualities as the Van Gogh but is specially designed for the inside bend of the stairs. Because the rail runs above the narrow side of the stairs, the broad section of the stairs remains completely free for other stair users.



Automatic folding hinge

Our automatic folding hinge is easily operated via the joystick and remote controls. Only 800mm of wall space is required from 1st nose to back wall.

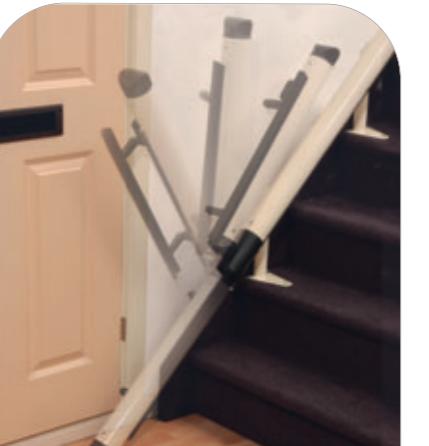

Handicare's tight-fitting drop nose rail can save having to fit a hinged track at the bottom of the stairs. A wall length of only 250mm is required.



Seat Features



Joystick

This easy to use control can be positioned on left or right arm

Key Switch

The key can be removed to prevent unauthorised use



Space Saving Design

The arm rests and foot support can be folded away easily, maximising the clear space on the stairs.



Safety Belt
Fitted as standard for your protection and security



Seat belt Extension
Available as optional extra



Manual Swivel
Push down lever allows seat to be swivelled enabling safe mounting and dismounting

Powered Swivel*
is available, operated by the joystick

Footrest Linkage*
Simultaneously folds up with seat



Off-centre Footrest*
A unique design available for tight winding steps
*optional

45° Downward Facing Seat

The Classic Seat is able to travel downward facing at 45°. This means that where a stairwell is particularly narrow, or a user has limited knee flexion, you can fit a Freelift Van Gogh rail with a 45° seat and special footplate. Minimum stair width 690mm (approx 27").

This could avoid the inconvenience of major building adaptations, or save on the cost of a through floor lift.


45° Travel

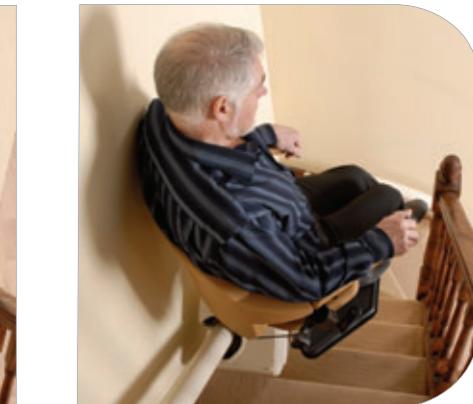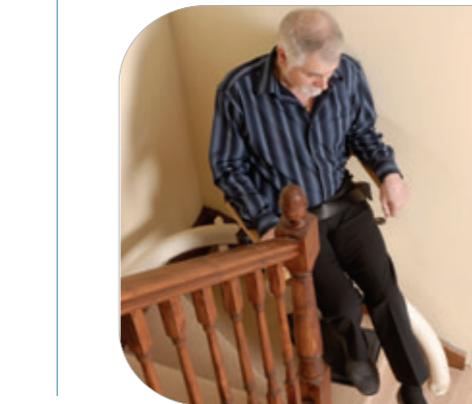


Rail end variations - top and bottom




90° top parking bend

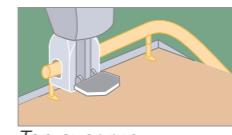


Top over-run

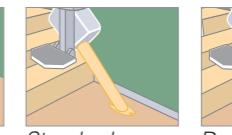

Standard

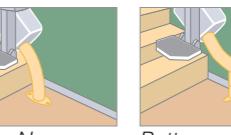


Drop Nose

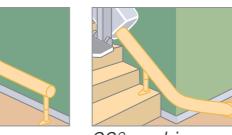

Bottom over-run



90° parking bend bottom

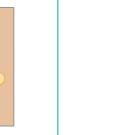

180° parking bend bottom


Installed in half a day without any construction work

We always attach great importance to installing your stairlift to the highest professional standards.

Experienced engineers complete the job in half a day, in most cases, without any alterations to your staircase.

Once a Handicare is installed the engineer will demonstrate and explain just how easy your stairlift is to use.

The arm rests and foot support can then be folded away, clearing as much space as possible for other stair users.

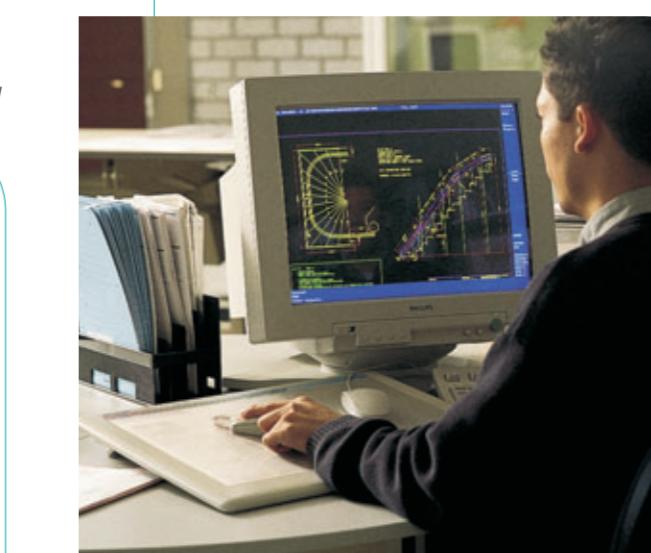

Why choose a Handicare Stairlift?

- You remain in your own familiar surroundings
- Safe, secure, easy transport up and down stairs
- Smooth and quiet operation
- Slim and stylish design for comfort and ease of use
- Chairs with stylish leather grain effect or comfortable fabric
- Installation within half a day and without damage to walls
- Full warranty and comprehensive after sales service