


Toilet seat frames and surrounds

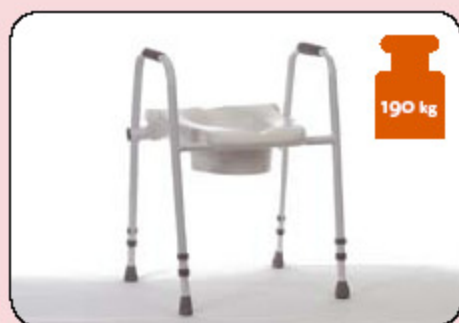

Merlin raised toilet seat frame

The adjustable height frame is made from powder coated tubular steel and features a contoured plastic seat which provides a high sitting position when placed over a toilet. The seat can be removed for cleaning. Cushioned handles provide extra support when transferring. Fitted with non-marking rubber ferrules to prevent slipping.

Merlin raised toilet seat frame 10380C

Technical specification

Product weight	5.2 kg	Seat height	460 - 610 mm (18- 24")
Width	555 mm (22")	Max user weight	190 kg (30 st)
Depth	495 mm (19.5")		



Aluminium toilet surround

This adjustable height toilet surround is constructed from anodized aluminium with non-corrosive fixings. Cushioned handles provide extra support when transferring. Fitted with non-marking rubber ferrules to prevent slipping.

Aluminium toilet surround 10397C

Technical specification

Product weight	1.9 kg	Arm height	620 - 695 mm (24.5 - 27")
Width	540 mm (21")	Max user weight	160 kg (25 st)
Depth	420 mm (16.5")		