

OPTIONS & ACCESSORIES

QUICKIE®





W

QUICKIE®

Real life has many faces. Life can be exciting, challenging, surprising. Whatever life has in store for you, it's how you perceive it that matters. Every single day offers unique chances and possibilities. The question is:

What's your **Attitude?**

Introducing the new Attitude range from Quickie; the Argon® / Argon® Ti / Neon® / Neon® Swing Away offer new possibilities and underline your uniqueness in both attitude and lifestyle without compromise.

This brochure outlines the different features and options the Attitude series offers. You will find an extensive range of accessories that complement this range, all of which can be installed on chair build or added to your existing Quickie chair by one of our Official Quickie Dealers.

Section 1 - Seating

Seat Slings	4
Seat Cushions	5

Section 2 - Backrests

Backrest Frames	6
Backrest Upholstery	8
Push Handles	9
Other Options	10

Section 3 - Armrests & Side guards

Armrests	12
Side guards	14

Section 4 - Legrests

Hangers	16
Footplates	18

Section 4 - Axle Plates

21

Section 5 - Castors & Rear Wheels

Castor Forks	23
Castors	24
Rear Wheels	26
Tyres	28
Handrims	30
Spoke Guards	31
Wheel Locks	34

Section 6 - Accessories

32

Please note the options displayed in this catalogue may not be available in your country. All options are subject to compatibility with other options you may have chosen. Always refer to the order form for option availability and compatibility.

Please note the reference codes in this brochure are NOT part numbers and cannot be ordered. Reference codes relate to the order form and are for reference purpose only.

Please refer to the relevant parts manual for part numbers.

SEATING

The extensive and versatile range of seating options from Quickie provides users with a seating solution to match their lifestyle so whatever Quickie chair they choose it is ready to meet their needs. For seating solutions for users with more demanding clinical requirements see our **JAY** Range.

Standard seat sling

The standard seat sling differs between the Argon, Argon Ti and the Neon. Both provide a rigid base for seat cushions and are available with a choice of utility bags located to the underside of the sling.

3002 - 1 Utility bag

3006 - 2 Utility bags

3005 - Catheter bag



Available Argon and Argon Ti



Available on Neon



Neon - Comes complete with a detachable padded cover in the same colour as the sling.

Lightweight seat sling

Made from a vented lightweight vented material, the lightweight seat sling provides a weight saving for active users where every gram counts.

Available with a choice of utility bags

3014 - 1 Utility bag

3017 - 2 Utility bags

3018 - Catheter bag



Available on the Argon Ti

3010 - Solid seat padded


Consisting of a wooden base and foam cushion this seat is available to suit all frame sizes.



Available on the Neon



SEAT CUSHIONS

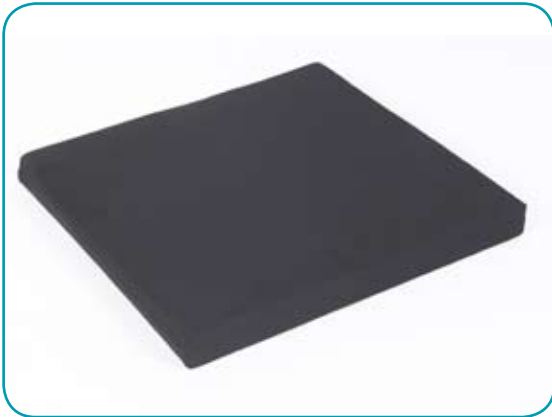

Available on Neon, Argon and Argon Ti

Seat cushion with black cover

The standard seat cushion, complete with removable washable cover, is available in foam composite or a softer latex material both available in a 3cm or 5cm cushion thickness.

Foam cushion - **3026**

Latex cushion - **3027**



Available on Neon, Argon and Argon Ti

Seat Cushion with water resistant / breathable cover

This seat cushion comes complete with a breathable inner cover that is water resistant whilst enabling air exchange when seated.

Like the standard cushion this comes in a 3cm or 5cm thickness with a choice of foam or latex cushion

Foam cushion - **3028**

Latex cushion - **3029**



Available on the Neon

Medical lamb skin cushion cover - **1404**

Ideally suited for users who require added warmth and a softer feel when seated, the medical lamb skin cover is completely removable and wraps around standard seat cushions.

Available in 3 sizes

35cm wide

40cm wide

45cm wide



BACKRESTS

The Quickie range of manual chairs offers a choice of backrests to suit any situation. Choose from a range of backrest frames and upholstery to suit users with mild to moderate clinical requirements. For seating and backrest solutions for users with more demanding clinical requirements see our extensive **JAY** Range.

BACKREST FRAMES

4017 - Fixed 90°, non angle-adjustable, non-folding backrest

Fixed at 90° this backrest is not adjustable and can not be folded. Ideally suited for users who want a rigid backrest with minimal play and movement when propelling.

Ideally suited for users looking to save weight and use this back frame in combination with a solid back



Available on Neon

4018 - Non-folding, angle-adjustable backrest

The backrest provides +8° (backwards) to -12° (forwards) of back angle adjustability in 5° increments.

This simple but effective design allows the backrest to be adjusted in minutes by relocating the fixing screws through a sequence of holes located at the rear of the frame.

Unique in design, this mechanism provides users with an adjustable backrest that is firm and rigid with no play.



Available on Neon, Argon and Argon Ti



Folding, angle-adjustable backrest - 4029

Using the same angle-adjustable mechanism as the non-folding, angle-adjustable back (4018) this back provides users with the same level of adjustment with the added benefit of allowing users to fold down the back in seconds using the quick release pull cord for ease of storage and transportation.



Available on Neon, Argon and Argon Ti



Double locking backrest - 4016

An addition to the folding backrest (4029), this allows the backrest to be locked whilst in the folded position with the flick of a lever.

This enables users to lift and move the chair ensuring the back stays in the folded down position.

Available on Argon and Argon Ti

UPHOLSTERY

Adjustable backrest upholstery

The Quickie range of backrest upholsteries all provide users with excellent support and flexibility through tension adjustable straps which can be adjusted and shaped to the user's back for greater comfort and better posture. Available in nylon (std), vented and leather.



4032 - Nylon / Standard - Available on Neon, Argon and Argon Ti

Available in a range of colours



4038 - Vented material - Available on Neon, Argon and Argon Ti

The vented material enables air exchange and heat dissipation for greater comfort when seated for long periods.



4045 - Leather - Available on Argon Ti

Made from nappa leather that is supple and soft to the touch but firm enough to maintain posture.



4025 - Solid backrest


Padded for extra comfort and suspended via the backrest tubes, this detachable solid backrest is ideally suited for users requiring extra support in the lumbar region.



Available on Neon



PUSH HANDLES

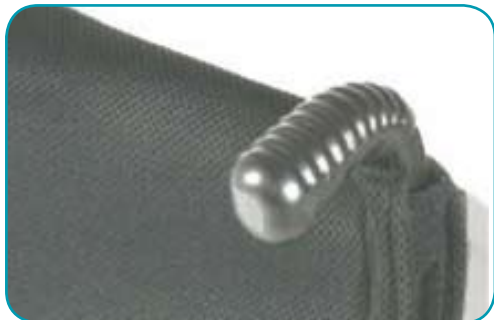

Available on Neon, Argon and Argon Ti

Short push handles - **4042**



Available on Neon, Argon and Argon Ti

Long push handles - **4022**



Available on Neon, Argon and Argon Ti

Fold down push handles - **4024**

Fold these push handles down with the click of a button for ease when storing or transporting the chair.



Available on Neon, Argon and Argon Ti

Height adjustable push handles - **4034**

Easy to adjust through a simple lever making life easy for users with different attendants and carers.

OTHER OPTIONS

4037 - Stabilising bar

Painted to match the frame colour and mounted onto the back tubes, the stabilising bar provides added rigidity to folding chairs. Ideally suited for use with tension adjustable backrest.



Available on Neon

4041 - Folding stabilising bar

Mounted onto the push handles the folding stabilising bar provides added rigidity to folding chairs, ideally suited for use with tension adjustable backrest.

Slide the sleeve from the centre and the folding stabilising bar folds with the chair.

This option is only available when the chair is fitted with height adjustable push handles.



Available on Neon

4027 - Headrest

The Quickie headrest provides both height and angle adjustability, and can be set asymmetrically to suit individual client needs.



Available on Neon



Available on Neon, Argon and Argon Ti

Jay hardware height adjustable push handles - 4043

This hardware kit includes the Quickie height adjustable push handles and contains everything that you need to fix any Jay® backrest to the push handles supplied.



Available on Neon, Argon and Argon Ti

Positioning belt - 4049

For added safety and security choose the Quickie positioning belt with its car style push to release button.



Available on Neon, Argon and Argon Ti

Seat depth growth kit - 4020

Inserted into the seat rail the seat depth growth kit provides users with an extra 5cm of seat depth often required when fitting third party backrests such as Jay® or for femoral growth.

This option cannot be chosen with the 90° fixed backrest (4017) and comes with a 70kg user weight limit



Available on Neon

Backrest cushion

Available in a choice of thickness and cushion material the Quickie backrest cushion provides extra comfort and support for users

3cm thick, foam cushion - 4035

5cm thick, latex cushion - 4036

ARMRESTS & SIDE GUARDS

The extensive range of armrests and side guards means whatever a user's clinical need or whatever their lifestyle demands we have an option to suit them

Fixed height armrest

The fixed height armrest comes complete with an integrated side guard painted to match the frame colour.

Using the trigger, users can flip the entire armrest back, rotate it to the rear of the backrest or completely detach it, aiding side transfers.

This armrest is available with (shown here) or without the armrest pad and T bar

4213 - Without pad

4214 - With pad



Available on Neon

Height adjustable armrest

Similar to the fixed height armrest this can be flipped back, swung around and completely detached.

Users can also adjust the armpad height using the push button located to the underside of the armpad.

This armrest is available with a choice of armpad sizes

4216 - Short pad (25cm)

4217 - Long pad (30cm)



Available on Neon



Padded Quickie-style armrest - 4265

The padded Quickie-style armrest is lightweight for active users and can be swung-away or completely detached to aid transfers.

Available on Neon, Argon and Argon Ti



Single post, height adjustable Quickie-style armrest - 4064

This detachable, height adjustable armrest comes complete with an integrated side guard and padded grip to help when transferring.

Available on Neon, Argon and Argon Ti



Quad release mechanism - 4246

For use with the fixed height (4013 & 4014) or height adjustable armrests (4016 & 4017) this adaption fixes directly onto the trigger of the armrest.

Available on Neon

SIDE GUARDS

4229 - Straight aluminium side guard

Aluminium and painted to match the frame colour, this straight side guard provides protection from debris, whilst improving comfort through a foam inlay, also aiding stability for the user.

Available as fixed (standard option) or with slide back option (4233)




Available on Neon, Argon and Argon Ti

4238 - Straight aluminium side guard with fender

Aluminium and painted to match the frame colour, this side guard extends and rolls over the wheel to provide added protection from debris, whilst improving comfort through a foam inlay, also aiding stability for the user.

Available as fixed (standard option) or with slide back option (4233), ideally suited for use with mountain bike wheels.

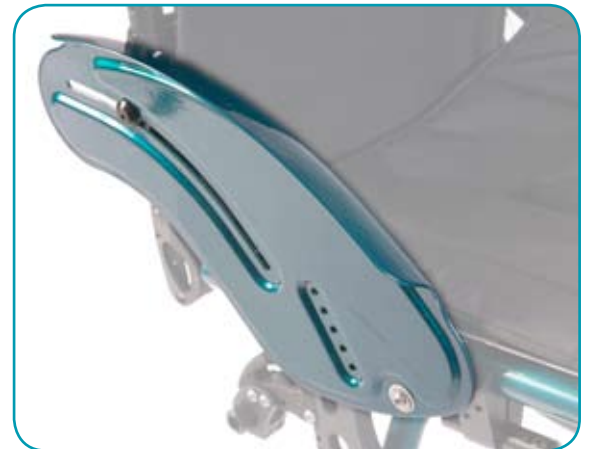

Available on Neon, Argon and Argon Ti

4232 - Straight carbon fibre style side guard with fender

This side guard provides added protection from debris, whilst improving comfort through a foam inlay, together aiding stability for the user.

Ideally suited for users looking to save weight.

Not available with slide back option



Available on Neon, Argon and Argon Ti



Small lightweight carbon fibre side guard - 4247

Made from lightweight carbon fibre this fixed side guard provides active users with protection and support without compromising on weight.

Only available on Argon Ti



Plastic detachable side guard - 4263

This black plastic side guard is lightweight and rigid providing users with support and protection and like all other side guards can be detached, aiding with side transfers.

Available on Neon, Argon and Argon Ti



Slide back option - 4233

Chosen in combination with the straight aluminium side guard (4229) or the straight aluminium side guard with fender (4238) this option allows the chosen side guard to be flipped back to aid with side transfers.

Available on Neon, Argon and Argon Ti

LEGRESTS

Fixed, swing-away or elevating, together with a wide range of footplates to suit every client preference.

HANGERS

Quickie-style hanger

Available in 70° or 90°, the Quickie-style hangers can be swung in or out to make room when transferring or completely removed for storage and transportation.

These hangers offer incremental lower leg length adjustment providing 22cm to 50cm lower leg length, dependent on the chair configuration.



Available on Neon Swing Away.

Sopur-style hanger

Available in 75°, the Sopur-style hangers can be swung in or out and completely removed. With new cambolt technology play is dramatically reduced.

With a larger release mechanism this hanger is easier to operate for people with reduced dexterity and strength.

Enables 22cm to 50cm lower leg length, dependent on the chair configuration.



Available on Neon Swing Away.




Available on Neon Swing Away

Elevating legrest

The Quickie elevating legrest can be easily adjusted by the flip of a lever.

Choose right hand (2129), left hand (2130) or a pair (2128) dependent on user requirements.

Supplied with the divided platform footrest, painted in the frame colour, this legrest can be swung in, out or completely detached.

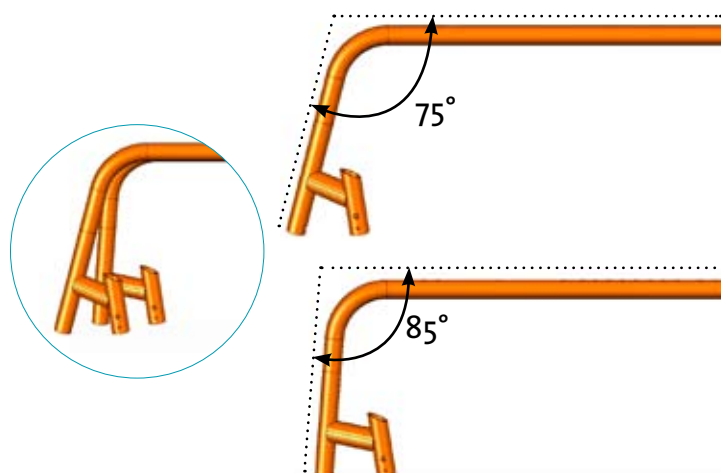

Fixed hangers / front frame


Two variations of front frame angle are offered on the Quickie Argon, Argon Ti and Neon Fixed Front.

75° - 1241

85° - 1270

Frame inset on fixed hangers

Available on Neon Fixed Front, Argon and Argon Ti, choose either a 0cm, 2cm or 4cm frame inset. The frame inset is the distance the front frame tapers inwards on each side providing users with a sleek stylish look, additional support around their legs and makes manoeuvring through tight spaces easier.

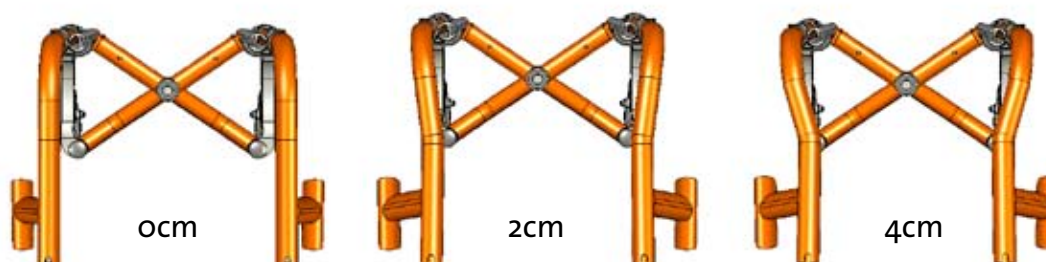

4cm not available on Neon Fixed Front

FOOTPLATES

2100 - Flip-up footplates

Our tried and tested Quickie flip-up footplates made from a composite plastic are both lightweight and rigid.

Used with swing-away hangers (Quickie-style shown here) users can flip the footplates up or swing the entire legrest out of the way.

Also available in an angle-adjustable version to suit more specific seating requirements - **2106**



Available on Neon Swing Away

2119 - Platform footrest, aluminium

Made from aluminium and painted to match the frame colour, this enhances the rigidity of any folding chair and will benefit users looking for a stable and rigid platform.

This footplate comes complete with a horizontal calf strap and has the ability to adjust the angle to suit individual user requirements and to be flipped up out of the way when required.



Available on Neon

2120 - Divided footplates

Similar to the platform footrest, the divided footplates are also constructed from aluminium and painted to match the frame colour.

Angle and height adjustable (via extension tubes) these footplates can also be flipped up for ease of transfer and swung away using the swing-away hangers.



Available on Neon Swing away, Argon



Tubular footrests - 2106

This lightweight aluminium footrest provides excellent chair rigidity and a solid platform.

Available on Argon and Argon Ti



Tubular footrests with flat cover - 2105

Like the tubular footrest with the addition of a snap-on flat cover.

Available on Argon and Argon Ti



Platform footrest, angle and depth adjustable - 2101

The platform footrest provides users with both angle and depth adjustability through the footplate whilst they are able to adjust the lower leg length through the extension tubes.

Available on Argon and Argon Ti

2131 - Platform footrest, angle and depth adjustable

Like the Argon platform footrest this provides users with both angle, depth and lower leg length adjustment.

This footrest can be flipped up for transfers or when the chair is folded.

Available in black composite or painted aluminium to match the frame colour.



Available on Neon

2160 - Auto-folding footboard

The auto-folding footboard provides users with a strong and rigid base that can be depth, height and angle adjusted to suit individual needs.

Available in black this footboard automatically folds up when the chair is folded making transportation and storage quick and easy.



Available on Neon Fixed Front

2133 - Foot positioning plates

The Quickie foot positioning plates provide users with added support, restricting movement.

For use on platform footrest.



Available on Neon Fixed Front , Argon and Argon Ti

High mount footrest option

The high mount footrest option is ideal for users with shorter lower leg lengths. Available in two styles; divided or platform. The footrest supplied differs depending on the chair.



Fine tuning and personal adjustments are obtained through a range of axle plates allowing users to adjust the centre of gravity and seat height to their own personal preference.

AXLE PLATES



Standard

Neon - Axle plates

The range of Neon axle plates has been designed to enable fast and easy adjustment of both seat height and centre of gravity, through a series of holes and bolts.

Each axle plate is unique in design and provides a different level of adjustment - For more information on the adjustability gained through each axle plate please refer to the order form.



Short



Amputee



Handbike



Suspension

Argon & Argon Ti Axle plates

The Argon and Argon Ti axle plates provide the user with fast and easy adjustment of both seat height and centre of gravity.

Unlike the Neon axle plates all of these axle plates provide infinite adjustment within a range - For more information on the adjustability gained through each axle plate please refer to the order form.



Argon - Standard



Argon Ti - Short & Standard



Argon & Argon Ti - Amputee



Argon - Handbike



Argon & Argon Ti - Suspension



Choose from an extensive range of castors, rear wheels and tyres from Quickie and other manufacturers, such as Schwalbe and Spinergy.

CASTORS
&
REAR WHEELS

CASTOR FORKS



Standard castor fork

Available on the Argon and Neon, the standard castor fork is available in a variety of sizes.

Optional quick-release castor axles are available.

Available on Argon and Neon



Argon Ti castor fork

Lightweight and stylish, the new Argon Ti castor fork is available in 98mm, 111mm and 123mm.

This castor fork is also available in an extra wide version to accommodate the new extra wide soft roll aluminium castor and is available with optional quick-release castors.

Available on Argon Ti

CASTORS

Pneumatic castor

For those users who prefer a softer ride over rough terrain, these air filled castors provide more comfort than the solid equivalent.

Available in 6" **5823** (shown here) and 7" **5838**.



Available on Neon

5812 - Solid castor

These solid polyurethane castors are maintenance free and have a very low rolling resistance for increased energy efficiency.

Available in 3", 4", 5", 6" and 7", smaller sizes offer a better turning radius and greater foot clearance.



Available on Neon, Argon and Argon Ti



Soft solid castor - 5823

These solid castors have a wider rim and are made from a softer composite providing a smoother ride whilst still been maintenance free.

Available in 5" and 6"

Available on Neon, Argon and Argon Ti



Aluminium rim castors - 5852

These low profile, maintenance free solid castors are made with a stylish aluminium rim and offer an extremely responsive ride .

Available in 4" and 5" and new extra wide version (5855)

Available on Neon, Argon and Argon Ti



Solid light-up castor- 5850

Stand out from the crowd with a set of light-up castors. These castors are maintenance free and are available in 3" and 4".

Available on Neon, Argon and Argon Ti

REAR WHEELS

5171 - Straight spoke wheel

These standard straight spoke wheels come in 20", 22" and 24", complete with quick-release axles.



Available on Neon, Argon and Argon Ti

5172 - Cross spoke wheel

These standard cross spoke wheels come in 20", 22", 24", 26" complete with quick release axles.



Available on Neon, Argon and Argon Ti



Design wheel - 5137

This stylish highly polished aluminium wheel has straight radial spokes and is available in 20" and 24"

Available on Neon, Argon and Argon Ti



Lightweight wheel - 5140

The lightest standard wheel on the market. These wheels provide 1.5kg weight saving over a pair of straight spoke wheels.

Not available with solid tyres (5168).

Available on Neon, Argon and Argon Ti



Everyday Spox - 5159

Made with a double-walled custom designed aluminium rim for strength, these wheels are extremely reliable and comfortable for everyday use.

Lightweight, reliable and maintenance free.

Available in 24", 25" and 26"

Available on Neon, Argon and Argon Ti

TYRES

5124 - Standard profile tyre - Pneumatic

This air filled tyre provides a cushioned ride with a stronger profile, providing better grip and traction than slick tyres.

Available in 22" and 24".



Available on Neon, Argon and Argon Ti

5125 - Standard slick tyre - Pneumatic

This air filled tyre offers a cushioned ride and a slick tread for low rolling resistance.

Available in 22", 24" (Schwalbe right run) and 26" (slick tyre)



26"



Available on Neon, Argon and Argon Ti

5118 - Lightweight Inova - Pneumatic

Air filled providing a cushioned ride and with a slick tread, these tyres are ideally suited for active users looking for a weight saving where every gram counts.

Available in 24" only with lightweight wheel.



Available on Neon, Argon and Argon Ti

5168 - Solid tyre

This black non-marking solid tyre has a smaller profile for low rolling resistance

Available in 24".

Not available with lightweight wheel.



Available on Neon, Argon and Argon Ti

5233 - 22" x 13 / 8" puncture proof - Solid

This black non-marking solid tyre with a stronger profile provides increased grip and traction.

Available in 22".

Not available with lightweight wheel.



Available on Neon, Argon and Argon Ti



Schwalbe Downtown - 5167

2grip



At each push and particularly when braking, the hand's constant contact with tyres can result in irritation and chafing. The Schwalbe Downtown with its patented 2 grip solution consists of completely smooth sidewalls, allowing hands to glide past the tyres without irritation.

Available in 24"

Available on Neon, Argon and Argon Ti

Schwalbe Speedrun - 5180



The new Speedrun folding tyre combines the requirements of handbikers, indoor sportsmen and everyday users. It is ultra light at just 190g and includes a highly effective puncture protection layer. The black rubber compound will not leave unsightly marks and so is also suitable for indoor use.

Available in combination with 25" Spinergy or Lightweight wheel

Available on Argon Ti

(((smartguard®)))



Schwalbe Marathon Plus Evolution - 5164

The newly developed puncture protection belt SmartGuard® – is extremely effective at preventing tube defects. The insert, made from highly elastic, special india rubber, offers durable resistance to penetration punctures. The Marathon Plus combines the comfort and easy rolling features of a pneumatic tyre with puncture resistance associated with solid tyres.

Available in 24"

Available on Neon, Argon and Argon Ti

2grip



Schwalbe Mountain Bike Wheel - 5162

The Mountain bike wheel option comes complete with a 24" mountain bike wheel, black powder coated handrim and Schwalbe 2 grip tyre.

This off-road tyre for wheelchairs is designed with sturdy lugs and large volume so it performs well off-road.

Available on Neon, Argon and Argon Ti



HANDRIMS



5144 - Aluminium silver powder coated handrim

Aluminium silver anodized handrim - **5142**

5163 - Vertical projection handrim

Coloured aluminium handrim - **5183**

5174 - Stainless steel handrim

Titanium handrim - **5122**

5121 - Titanium clear powder coated handrim

Supergrip handrim - **5120**



Add a personal touch with our range of funky spoke guards to suit 22" and 24" rear wheels. Transparent spoke guards are also available (not shown here) to suit 22", 24" and 26" rear wheels.

SPOKE GUARDS



5135 - Cheese



5129 - Bibbo



5134 - Galactica



5127 - Rotation



5130 - Fantasy



5128 - Graphic

WHEEL LOCKS

6009 - Standard wheel lock

Spring-loaded for ease of use, when engaged the lever lies out of the way making side transfers easy.



Available on Neon, Argon and Argon Ti

6008 - Knee lever wheel lock

Spring-loaded for ease of use and mounted at knee level for easy access, this wheel lock is suitable for most users.



Available on Neon, Argon and Argon Ti

6016 - Extension handle

This extension handle makes operating the knee lever wheel lock easier for users with restricted movement.

Available in combination with knee lever brake



Available on Neon, Argon and Argon Ti



Available on Neon, Argon and Argon Ti

Compact wheel lock - 6003

Ideally suited for users with reduced dexterity .



Available on Neon, Argon and Argon Ti

Grade aid - 6005

Designed to prevent roll-back on an incline, the grade aid provides security and peace of mind when travelling up steep inclines.



Available on Neon, Argon and Argon Ti

Scissor wheel lock - 6004

The spring-loaded scissor wheel lock can be engaged with one simple movement.

ACCESSORIES

1205 - Crutch holder

The easy to install crutch holder keeps crutches out of the way but still in reach.

Available in combination with tip assist (1230)



Available on Neon, Argon and Argon Ti

1244 - Anti tip tubes

Prevent backwards tipping. Can be foot operated by the attendant.



Available on Neon, Argon and Argon Ti

1245 - Quickie-style anti-tip tubes



Available on Neon, Argon and Argon Ti



Available on Neon, Argon and Argon Ti

Transit wheels - 1247

Specifically designed to allow users to access tight spaces by removing the rear wheels.



Available on Neon, Argon and Argon Ti

Tip assist - 1230

This option allows attendants to easily tip the chair back to mount kerbs and steps.

Available in right hand, left hand or as a pair.



Available on Neon, Argon and Argon Ti

Frame protection - 1281

Protect the frame from everyday dents and chips with this easily applied protection sleeve. Fits onto the front frame or hangers.

1403 - Air pump

Bike style air pump.



Available on Neon, Argon and Argon Ti

1402 - High pressure air pump

High pressure, quality track pump.



Available on Neon, Argon and Argon Ti

4006 - Quickie back pack

The stylish Quickie back pack hooks directly onto the push handles and is ideal for storage.



Available on Neon, Argon and Argon Ti



Mobile phone holder - 2113



Bag caddy - 2118

Can be used for carrying laptops or briefcases.



Net bag - 1405

Available on Neon, Argon and Argon Ti



Sunrise Medical Limited
High Street
Wollaston, Stourbridge
West Midlands, DY8 4PS
T +44 (0) 1384 446622
F +44 (0) 1384 446644
www.sunrisemedical.com

INNOVATE.™
BUILD FOR THE FUTURE