



Salsa M



Salsa M

Excellent indoor manoeuvrability with great outdoor performance



Excellent indoor manoeuvrability


- Quickie mid wheel drive technology delivers an ultra-compact turning radius for superior indoor manoeuvrability
- Narrow width of 60cm enables easy access through doors and confined spaces
- Very low seat height (from 40cm, and from only 43cm with powered tilt) makes access under tables and into vehicles easy

Great outdoor performance

- Patented Quickie all-wheel independent suspension ensures a smooth and stable ride over rough or uneven terrain
- Quickie anti-pitch technology eliminates pitching on slopes, ensuring your safety at all times
- Effortlessly climbs up to a 70mm kerb
- 50Ah gel batteries and powerful high-torque motors provides up to 26km (16 miles) of driving range at speeds of 6 or 10kph (4 or 6mph)

Fully adjustable seating system

The Salsa M seating system with width and depth adjustment ensures a perfect fit today and tomorrow.

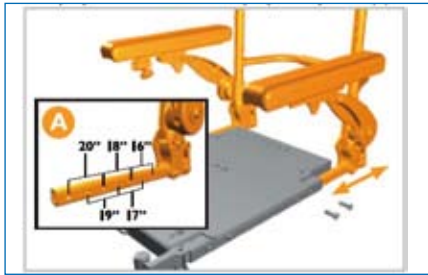

Flip-up armrests and foot rests for easy accessibility

- Seat width and depth adjustment from 41 to 51cm
- Backrest height adjustable in three positions from 45.5 to 50.5cm
- Five adjustable seat heights and four seat angles
- Easy adaptation for alternative seats, backs and cushions (e.g. standard back sling, Jay, third party seating)



NEW: Low push force, emergency freewheel allows the Salsa to be moved very easily when it is in freewheel



compatible

Lift-off backrest for easy transportation and storage



Full range of powered options available

Available powered options:

- C.O.G tilt-in-space (0° to 30°)
- Lift and tilt (30cm lift, 0-30° tilt)
- Recline
- ELR/ALR powered leg rests
- Easy adaptation of powered options if required (upgrade kit)

Optional retro-fit kits:

A full range of upgrade kits can be added to your chair at a later date if your needs change

Advanced electronics

VR-2 electronics as standard

- 70 Ah controller as standard. 90 Ah heavy duty version optional
- Dual control (option)

R-net expandable controls (optional)

- Provides simple programming options
- Shortcut to your favourite functions at the touch of a button thanks to Quickie's patented assignable button technology
- Control up to 6 powered functions through your joystick
- Drive profiles can be tailored to your personal requirements
- Control your PC with your joystick via Bluetooth® (option)
- Dual control (option)

Easy to service



Salsa M

Freedom indoor and out



Technical Data: Seating

Seat angle:	Adj. 0°, 3°, 6°, 9°	
Seat depth:	Adjustable from 41 to 51cm	16 - 20"
Seat width:	Adjustable from 41 to 51cm	16 - 20"
Seat height:	40 - 45cm (no tilt), 43 - 48cm (powered tilt), 48 - 53cm (lift & tilt)	15.7 - 17.7" (no tilt), 17 - 19" (powered tilt), 18.5 - 20.5" (lift & tilt)
Backrest angle:	Adj. -4°, 0°, 4°, 8°, 12°, 16°	
Backrest height:	Adj. 46, 48 and 51cm	18, 19 and 20"
Backrest width:	Adjustable from 41 to 51cm	16 - 20"
Lower leg length:	33 - 51cm	13 - 20"
Legrest angle:	70°, 90°, ALR, ELR	
Armrest height:	Adj. from 22 to 31cm	8.7 to 12"
Armrest depth:	Adj. from 31 to 41cm	12 to 16"

Specification

Overall width:	60cm	23.6"
Overall length: (including legrests)	110cm	43"
Chair weight*:	113 - 165kg	249 - 363lbs
Turning circle:	120cm	47"
Battery capacity	2 x 50 Amph	
Range:	26 km (ISO 7176-4)	16 miles
Speed:	6 kph and 10 kph	
Kerb climbing ability:	7cm	2.7"
Maximum user weight:	140kg	22 stone
Max safe slope:	8° (14%)	
Ground clearance:	5cm	2"
Transit tested to ISO 7176 Part 19 and passed:	Yes	

* Please note: Total weight of the powerchair depends on the configuration
113kg (base chair) to 165kg (base chair plus solid tyres, power lift, tilt, recline, headrest, comfort seating, R-net controls and lights)



Sunrise Medical Limited
High Street, Wollaston
Stourbridge
West Midlands
DY8 4PS

Tel: +44 (0)1384 44 66 66
Fax: +44 (0)1384 44 66 44
www.SunriseMedical.com