

QUICKIE®

Power



Salsa®



**Compact yet powerful...
...at home outside as well as inside!**

Configured Salsa® with Comfort Seating and powered elevating legrests



QUICKIE®

Compact dimensions for INDOOR mobility...



Low seat height for easy table access

- with 12" wheels at 0° manual tilt = 14.5" (37 cm)
- with 14" wheels = 15" (38 cm)
- with 30° power tilt module = 16" (41 cm)
- 4 x 0.8" (2 cm) seat height adjustments



Narrow chair width easily accesses narrow doors and confined spaces

- with 12" wheels chair width = 23" (58 cm)*
- with 14" wheels chair width = 24" (61 cm)*
- Flip-back lighting (option) does not impact width



Unique 11" (28.5 cm) base unit height with 60 Amp batteries as standard

- New battery, wiring and electronic layout for compact base
- Adjustable suspension for optimum comfort and handling
- Suspended frame for improved handling and traction control

... and performance for OUTDOOR freedom



3" (8cm) ground clearance for thresholds, ramps, kerbs and rough terrain with 14" wheels


- 14" wheel version climbs 4" (10cm) kerbs with kerb climber (std)
- 12" wheel version climbs 3" (8.5cm) kerbs with kerb climber (std)
- Easily mounts 9" (23 cm) dual ramps



Compact 60 Amp gel batteries and electronics

- Provides range of up to 20 miles (32km) [ISO 7176 Part 4]
- Maximum user weight 22 stone (140 kg)
- Safe to transport in vehicles, trains and planes
- Fully crash test approved (including J3 back and options)

* with larger seat sizes overall chair width may be increased

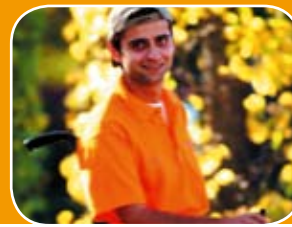

Innovative compact base = excellent **INDOOR** mobility + **OUTDOOR** freedom



Salsa 12" rear wheel option = 23" (58cm) chair width



Standard NHS PACK 3 Salsa® – 14" rear wheels and kerb climber (not shown above) are standard



Fully adjustable seating system with retro-fittable options for easy assessment, refurbishment and servicing



ONE fully adjustable seat with more than TWO MILLION seating combinations as standard!

- Only four tools required for adjusting complete seating system
- Five seat widths and five seat depths
- Five seat heights and four seat angles
- Three hanger depths plus adjustable footplates
- Full height, depth and angle adjustment on flip-back armrests
- Six backrest angles and three backrest heights – complete with adjustable push handles and width and tension adjustable upholstery



Easy transfer

- Flip-back armrests and controller
- Detachable, swing-away legrests
- Low seat to floor height
- Easily removable skirt guards




Easy to service

- Pull out battery tray for servicing with client in chair
- Easily removable shroud with accessible wiring/electronics
- Only four bolts to separate seat system from base



Safety on slopes, ramps and in transport

- Fully transit tested (ISO 7176 Part 19) with six point tie downs
- Quick release, removable backrest for easy transport
- Centre of Gravity (COG) seat adjustment for optimised stability
- Central batteries to improve stability
- Tested for use with 9" (23 cm) dual ramps

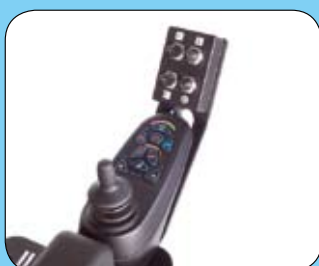

Easy to adapt retro-fit kits; ideal for refurbishment / upgrade

- Full range of retro-fittable options and accessories
- Up to four actuators and lighting (options)
- Compatible with new J3 backrests and positioning options
- Options can be ordered and delivered as standard kits

* For option and kit availability please see Salsa® configured and kit order forms





Full range of options available OR separately as retro-fittable option kits!



- Quickie Salsa® is easy to adjust, service, upgrade and downgrade to changing customer requirements
- Most option choices for the configured Salsa® are ALSO available as retro-fit kits for upgrade. Either order your Salsa® fully configured for specific clients or stock the standard Salsa® alongside retro-fittable kits for operational flexibility.
- Each retro-fit kit is delivered packaged and labelled with all the components required to fit an option, including the fitting instructions.
- See the Salsa® kit order form for further details and availability

Quickie Salsa®	KIT CODE
Upgrade kit: pair of 14" solid drive wheels + pair of 9" solid castors	060552
Upgrade kit: pair of 12" pneumatic drive wheels + pair of 8" solid castors*	060539
220V 8 Amp automatic charger with UK power lead	060540
Four actuator control*	060541
Dual attendant control kit*	060542
Manual backrest recline kit 16"-18" (41 - 46 cm) seat depth*	060544
Manual backrest recline kit 19"-20" (48 - 51 cm) seat depth*	060543
Power seat tilt*	060545
Power backrest recline 16"-18" (41 - 46 cm) seat depth*	060546
Power backrest recline 19"-20" (48 - 51 cm) seat depth*	060547
Power seat lift & tilt*	060548
Standard 2" cushion for 16"-18" (41 - 46 cm) width	060246-001
Standard 2" cushion for 19"-20" (48 - 51 cm) width	060246-002
Standard 70° legrests with adjustable footplates (pair)	060511
Manual elevating legrests with angle adjustable footplates (pair)	060512
Power elevating legrests*	060538
Comfort Seating Kit 16" (41 cm) width (seat & back and fixing kit)*	060549
Comfort Seating Kit 18" (46 cm) width (seat & back and fixing kit)*	060550
Comfort Seating Kit 20" (51 cm) width (seat & back and fixing kit)*	060551
J3 backs and positioning options	see J3 order form
JAY cushions	see JAY order form
Kerb climber [climbs 4" (10 cm) kerbs with 14" rear wheel]	060260
Standard headrest kit	060504
Crutch holder kit	060334
Swing away tray kit (LH or RH)	060494
Built 4 Me Specials	Call Tel: 01384 446521

* For option and kit availability please see Salsa® configured and kit order forms

Standard Salsa® seating system; ONE seat = over two million seating combinations!

inches / cms		Seat width				
		16 / 41	17 / 43	18 / 46	19 / 48	20 / 51
Seat depth	16 / 41					
	17 / 43					
	18 / 46					
	19 / 48					
	20 / 51					

4 x Seating angles

Seating angles	0°	3°	6°	9°
degrees				

5 x Seat heights

12" rear wheel	14.5 / 37	15 / 38	16 / 41	17 / 43	17.5 / 45
inches / cms					
14" rear wheel	15 / 38	15.5 / 40	16.5 / 41	17.5 / 44	18 / 46

5 x Centre of gravity settings

COG	-2 / 5	-1 / 2.5	0	+1 / 2.5	+2 / 5
inches / cms					

3 x Hanger depths

Hanger depths	+1 / 2.5	+2 / 5	+3 / 7.5
inches / cms			

3 x Backrest heights

Backrest height	18 / 46	19 / 48	20 / 51
inches / cms			

6 x Backrest angles


Backrest angles	-4°	0°	4°	8°	12°	16°
degrees						

Infinite armrest depth settings

Armrest depth	> 12 / 31	< 16 / 41
inches / cms		

Infinite armrest height settings

Armrest height	> 9 / 22	< 12 / 31
inches / cms		

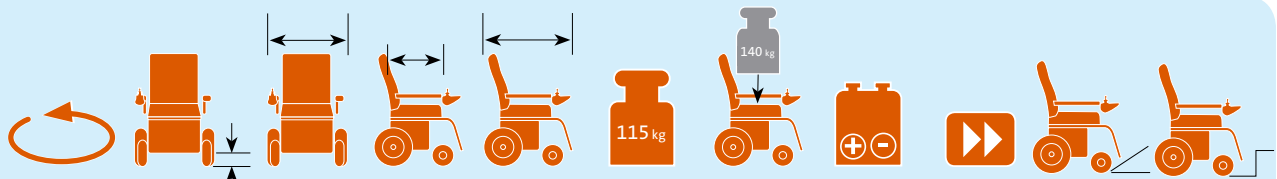

	Turning*	Ground Clearance	Chair Width	Chair Length no Legrests	Chair Length std Legrests	Chair Weight**	Max User Weight	Batteries	Speed	Max Angle of Climb***	Kerb Climb****
12" wheels	33.5" / 47" / 85 / 120 cm	2.8" / 7 cm	from 23" or 58 cm	31.5" / 80 cm	45" / 114 cm	114.5 kg	22 stone 140 kg	60 amp	4 mph 6 kph	25° / 10.6° / 8°	3.3" / 1.6" / 8.5 / 4 cm
14" wheels	33.5" / 47" / 91 / 120 cm	3" / 8 cm	from 24" or 61 cm	31.5" / 80 cm	45" / 114 cm	114.5 kg	22 stone 140 kg	60 amp	4 mph 6 kph	25° / 10.6° / 8°	4" / 2" / 10 / 5 cm

* Turning radius / Turning space for 3 point turn ** chair weight with batteries included. Salsa® weighs (77.5 kg) without batteries. Total chair weight is dependent on options chosen
 *** without client in chair / static with client in chair [ISO7176 Pt1] / dynamic with client in chair [ISO7176 Pt2] **** with kerb climber / without kerb climber

STLPS™
 Sunrise Training and Education Programmes
KNOWLEDGE FOR THE FUTURE

SUNRISE MEDICAL

Sunrise Medical Limited
 High Street
 Wollaston, Stourbridge
 West Midlands DY8 4PS
 T +44 (0) 1384 446666
 F +44 (0) 1384 446644
 www.sunrisemedical.com