

M A X I C O M F O R T TM S E R I E S


Rise



New!

MAXI
COMFORTTM
Ultimate Recline Technology


Softouch
with STAIN DEFENSE

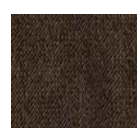

Hazelnut



Sterling

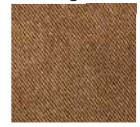

Copper

MaxiComfort Cloud

Whether you're looking for a lift chair to spend most of your day and night in, or you're looking for a chair that will rejuvenate you after a long day at work, look no further than the MaxiComfort Cloud – the most comfortable lift chair we've ever built! With an overstuffed, oversized biscuit style back, the Cloud features a unique seat and chaise pad design that cradles you in comfort.



Rise.

Roll.

Rest.

Ride.

GOLDEN
technologies

YOUR LIFE IN MOTION

MaxiComfort Cloud

M A X I C O M F O R T ™ S E R I E S

Model Number	PR-510
Size	Medium
Weight Capacity	375 LB
Recline Positions	MaxiComfort
Standard Fabrics	3
Fabric Options	5
Chaise Seat	Yes
Back Type	Biscuit
Head Rest Cover	No
Arm Covers	No
Battery Back Up	Yes
Overall Width	38"
Overall Height	45"
Width Between Arms	20"
Floor to Top of Seat	20"
Seat Depth	21.5"
Seat to Top of Back	29"
Seat to Top of Arm	7"
Floor to Back Edge of Seat Raised	28"
Floor to Front Edge of Seat Raised	24"
Seat Angle Raised	15°
Floor to Top of Foot Rest Reclined	24.5"
Floor to Top of Back Reclined	24"
Back Angle Closed	100°
Back Angle Reclined	170°
Foot Rest Angle Closed	90°
Foot Rest Angle Opened	180°
Distance Required from Wall Reclined	26"
Extension Reclined from Back of Seat to End of Foot Rest	44"
Life Time Warranty, Chair Frame, Lift Frame, Recline Mechanism	Yes
Warranty on Electrical parts	3 YR
Pro Rated Years 4 Thru 7	Yes
Warranty on Mechanical Labor	3 YR
Warranty on Electrical Labor	1 YR
Removable Back Option Available	No
Smart Tek System Standard	Yes
Customizable *	No
Left Hand Control Option Available	Yes
Extra Pocket Option Available	Yes
Heat & Massage Option Available	Yes
Foot Rest Extension Option Available	Yes
Straight Lift Option Available	No

* See Customize Chart for Dimension limits
 Sizes may vary depending on fabric, filling material, upholstery and carpet thickness. All measurements taken on concrete floor using levels and metal rulers.

Hand-blown, Overstuffed, Conjugated Microfiber Dacron Fill
SUPER SUPPORTIVE

Northern Maple Wood
 • Double-doweled at Stress Points
 • Reinforced Miter Blocks
 • Stapled and Glued
STRONGEST IN THE INDUSTRY

2" of Block Foam
EXTRA COMFORT

Quick Release Hand Control
EASY DISCONNECT

Rear Landing Legs
EXTRA STABILITY

5 Heavy Gauge Springs
LONG LASTING FOR INCREDIBLE COMFORT

Widest Fabric Selection in Industry
FITS IN ANY HOME DÉCOR

10 Gauge Steel Scissor Mechanism
STRONG & LONG LASTING

Adjustable Anti-Skid/Scrape Feet
STABILITY & SUPPORT



501 M shown



501 M shown



501 M shown



501 M shown

Golden Craftsmanship

Each luxury lift & recline chair is hand crafted in our facility using decades of motion furniture experience.

Added Comfort and Safety

Designed to perform year after year using advanced technology backed by the best warranty.

Smart Engineering

Smooth, quiet lifting system offers a secure lift for the best performance.



Smart Tek

Our luxurious lift and recline chairs are equipped with electronics to help assure a secure lift every time you need it.



Battery Backup System

Batteries not included. Battery Backup System is designed to bring the chair from the reclined to the seated position in the event of a power outage.



Literature is current at time of printing. Golden Technologies reserves the right to make changes to the product or literature at any time.



Complies with UL 73 and CSA C22.2 No. 68

Listing # E112417

